

HLT158 User Manual
Material 80G woodfree paper
8 Page Saddle Stitch Booklet
Two Signatures
Process: 1c/1c
Folded size: 125mm H x 125mm W




NIOSH N95 Particulate Respirator

USER INSTRUCTIONS

FT-N010
FPK1267rB

Thank you for your purchase!

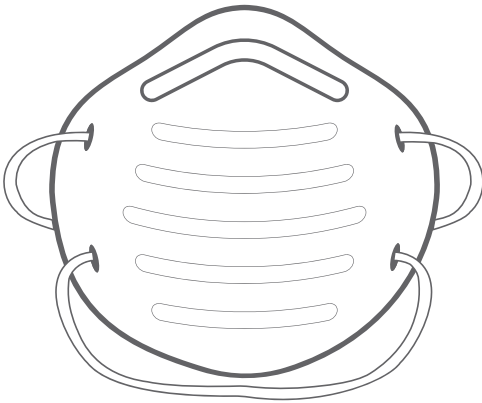

FLTR, Inc
www.FLTR.com
Toll free: 1-888-525-FLTR (3587)
(Monday - Friday 9am - 4:30pm PST)


WARNING:
This respirator helps protect against certain particles. **Misuse may result in sickness or death.** Materials which may come into contact with the wearer's skin could cause allergic reactions to susceptible individuals. For proper use, consult these Use Instructions. For questions or technical assistance, call FLTR toll free at 1-888-525-FLTR (3587).

- USE LIMITATIONS:**
- Do **NOT** use the respirator to enter or stay in a contaminated area under the following circumstances:
 - Atmosphere contains less than 19.5% oxygen
 - Atmosphere contains oil aerosols
 - For protection against gases or vapors
 - When contaminants or contaminant concentrations are unknown or immediately dangerous to life or health
 - When concentrations or contaminants exceed maximum use concentrations in applicable OSHA standards, applicable government regulations, or when they exceed 10 times the PEL (Permissible Exposure Limit)
 - For sandblasting, paint spray operations, asbestos
 - Do **NOT** modify or misuse the respirator.
 - Do **NOT** use the respirator if you have a beard or other facial hair that interferes with direct contact between the face and the edge of the respirator, or any other conditions that may prevent a good face seal.

SPECIAL INSTRUCTIONS

MODEL FT-N010 - N95 Filtering Facepiece Respirator has been manufactured by Suzhou Fangtian Industries Co., LTD. for FLTR, Inc. under TC-84A-7701

- FITTING INSTRUCTIONS:**
- Hold the respirator in your hand with the adjustable nose piece at your fingertips and let the straps hang downwards below your hand (fig. 1).
 - Place the respirator on your face, covering your nose and mouth completely. Pull both straps over your head. The upper strap should be placed on the top back of your head, above the ears. The lower strap should be placed on the back of your neck, below your ears (fig. 2).
 - Use two hands, bend the adjustable nose piece to conform to the shape of your nose (fig. 3).
 - A fit check should be performed each time a respirator is worn.
 - Cover the respirator with both hands and exhale quickly (fig. 4). If there is air leakage around the adjustable nose piece, repeat Step 3 to better conform the nose piece to the shape of your nose.
 - If air leakage is detected anywhere along the edge of the respirator, adjust the straps to achieve a better fit.
 - If air leakage cannot be eliminated, do not use the respirator and do not enter a contaminated environment.

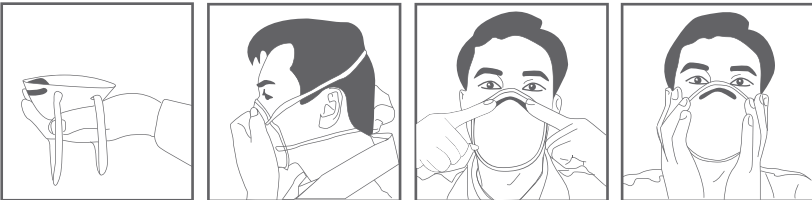

(Fig. 1) (Fig. 2) (Fig. 3) (Fig. 4)

USE INSTRUCTIONS:

- Failure to follow all instructions and limitations concerning the use of the respirator and/or failure to wear the respirator for the entire duration of exposure to contaminants can seriously reduce the performance characteristics of the respirator and lead to illness, injury, or death.
- Before use, the wearer must be educated in the correct use of the respirator in accordance with applicable safety and health standards.
- OSHA standards including 29 CFR 1910.134(e) (5) and CSA standard Z94.4-93 require that the wearer be fit-tested. Before occupational use of this respirator, a written respiratory protection program must be implemented meeting all local government requirements. In the United States, employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing.
- Discard the respirator and replace with a new one if:
 - Excessive clogging of the respirator causes breathing difficulty
 - The respirator becomes damaged
- Leave the contaminated area if dizziness, irritation, or other distress occurs.

DAMAGE INSPECTION

Check for holes in breathing zone. Leave contaminated area immediately if respirator is damaged, or if any physical distress such as breathing difficulty, dizziness or nausea occurs.



FLTR, Inc.
2821 Northup Way, Suite 200
Bellevue, WA 98004 USA
Telephone: 1-888-525-FLTR (3587)



THIS RESPIRATOR IS APPROVED ONLY IN THE FOLLOWING CONFIGURATION:

TC-	PROTECTION ¹	RESPIRATOR	CAUTIONS AND LIMITATIONS ²
		FT-N010	
84A-7701	N95	X	ABCJMNOPS

1. Protection

N95 - Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil; time use restrictions may apply.

2. Cautions and Limitations

- Not for use in atmospheres containing less than 19.5 percent oxygen.
- Not for use in atmospheres immediately dangerous to life or health.
- Do not exceed maximum use concentrations established by regulatory standards.
- Failure to properly use and maintain this product could result in injury or death.
- All approved respirators shall be selected, fitted, used and maintained in accordance with MSHA, OSHA and other applicable regulations.
- Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
- Refer to User's Instructions, and / or maintenance manuals for information on use and maintenance of these respirators.
- NIOSH does not evaluate respirators for use as surgical masks.
- Special or critical User's Instructions and/or specific use limitations apply. Refer to User's Instructions before donning.

NIOSH APPROVED: N95 PARTICULATE RESPIRATOR

95% EFFICIENT AGAINST PARTICULATE AEROSOLS FREE OF OIL.
THIS RESPIRATOR HELPS PROTECT AGAINST CERTAIN PARTICULATE CONTAMINANTS BUT DOES NOT ELIMINATE EXPOSURE TO OR THE RISK OF CONTRACTING DISEASE OR INFECTION. MISUSE MAY RESULT IN SICKNESS OR DEATH.

HLT158



NIOSH N95 Particulate Respirator

USER INSTRUCTIONS