

130mm

200mm

140mm

130mm

140mm

130mm

HENOX® 109013F




Particulate Respirator

20

NIOSH Approved Product

LOT NO.



**HENOX® SAFETY GROUP**  
NO 57 1ST + 2ND FLR, JALAN USJ 1/31  
TAMAN SUBANG PERMAL  
47600 SUBANG, JAYA, SELANGOR D.E. MALAYSIA  
Tel: 00603-8024-2906  
[www.henox-mkt.com](http://www.henox-mkt.com)

This respirator is only approved in the following configuration:			
TC	Protection <sup>1</sup>	Respirator	Cautions and Limitations <sup>2</sup>
84A-7867	N95	109013F	ABCJMNOPS

1. Protection

N95 - Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil. Time use restrictions may apply.

2. Cautions and Limitations

A - Not for use in atmospheres containing less than 19.5% oxygen.  
B - Not for use in atmospheres immediately dangerous to life or health.  
C - Do Not exceed maximum use concentrations established by regulatory standards  
J - Failure to properly use and maintain this product could result in injury or death.  
M - All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA and other applicable regulations  
N - Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.  
O - Refer to users instructions, and/or maintenance manuals for information on use and maintenance of these respirators.  
P - NIOSH does not evaluate respirators for use as surgical masks.  
S - Special or critical User's Instructions and/or specific use limitations apply. Refer to User's Instructions before donning.

S - Special or Critical User's Instructions: The model 109013F filtering facepiece respirator has been manufactured by Approval Holder Makrite for HENOX MARKETING SDN BHD under TC-84A-7867.

**SERVICE TIME LIMITATION**

1. Do not use the same respirator for more than 8 hours, continuous or intermittent, in dirty workplaces that could result in high filter loading.  
2. Leave contaminated area and remove respirator if breathing becomes difficult or proper fit cannot be obtained or the respirator is damaged, soiled, or distorted. Do not use the same respirator again once removed.  
3. The maximum service time of this respirator is five (5) consecutive calendar days (including days of non-use), beginning from the first day of use.

The model 109013F filtering facepiece respirator has been manufactured by Approval Holder Makrite for Company Henox under TC-84A-3323.

HENOX® 109013F



Particulate Respirator

LIMITED WARRANTY  
IMPORTANT NOTICE TO PURCHASER:


This limited warranty is made in lieu of the warranties of merchantability, fitness for particular purposes and all other warranties, express or implied. There are no other warranties which extend beyond the description on the face hereof. The physical standards and specifications of Henox will be met by products sold.

Exclusive Remedies: damages for the breach of this limited warranty are limited to the replacement of such quantity of Henox, products proved to be defectively manufactured. Except as provided above, Henox shall not be liable or responsible for any loss, damage, or liability, direct, indirect, incidental, special, or consequential, arising out of sale, use, or misuse, or the inability to use products by the user.

WARNINGS:

1. Weaver must read fitting instructions and use limitations before use.
2. This respirator is not guaranteed to protect against the risk of contracting a disease or infection.
3. Before occupational use of this respirator a written respiratory protection program must be implemented meeting all the local government requirements. In the United States employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing.

1

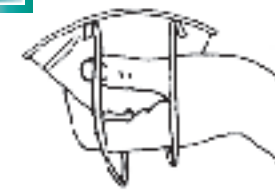

FITTING INSTRUCTIONS:

**Step 1**  
Hold respirator in your hand as shown, with the nosepiece at your fingertips and the headbands hanging freely below your hand.

2



**Step 2**  
Place respirator under your chin with the nosepiece up. Pull the top headband to top back of head.

3



**Step 3**  
Pull the bottom headband overhead and place around neck below ears.

4



**Step 4**  
Adjust headbands for comfortable fit. Push the nosepiece inward with your fingertips to conform to the shape of your nose.

**Caution:**  
The respirator should be fit checked for each wearing. Place both hands completely over the respirator and exhale. If air leaks around nose or the edges, adjust the nosepiece and/or headbands until a good fit is achieved.

