

150 MM

130 MM

150 MM

130 MM

175 MM

HF0904




20 Pieces N95 Disposable Respirator

HF0904




20 PIECES N95 Disposable Respirator





Apex Tool Group, LLC, 910 Ridgeway Road,
Suite 200, Sparks, MD 21152, USA
833-728-2866
www.satatools.com



This respirator is approved only in the following configuration:

| TC- | Protection¹ | Respirator | Cautions and Limitations² |
|----------|-------------|-------------|---------------------------|
| 84A-7867 | N95 | HF0904 X | ABC MNOPS |

1. Protection

N95-Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil; time use restrictions may apply

2. Cautions and Limitations

A - Not for use in atmospheres containing less than 19.5% oxygen
B - Not for use in atmospheres immediately dangerous to life or health
C - Do not exceed maximum use concentrations established by regulatory standards
J - Failure to properly use and maintain this product could result in injury or death
M- All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA and other applicable regulations
N - Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer
O - Refer to user's instructions, and/or maintenance manuals for information on use and maintenance of these respirators
P - NIOSH does not evaluate respirators for use as surgical masks
S - Special or critical User's Instructions and/or specific use limitations apply. Refer to User's Instructions before donning.

S - Special or Critical User's Instructions: The model HF0904 - filtering facepiece respirator has been manufactured by Approval Holder Makrite for Apex Tool Group, LLC, under TC-84A-7867


SERVICE TIME LIMITATION

1. Do not use the same respirator for more than 8 hours, continuous or intermittent, in dirty workplaces that could result in high filter loading.

2. Leave contaminated area and remove respirator if breathing becomes difficult or proper fit cannot be obtained or the respirator is damaged, soiled, or distorted. Do not use the same respirator again once removed.

3. The maximum service time of this respirator is five (5) consecutive calendar days (including days of non-use), beginning from the first day of use.

HF0904

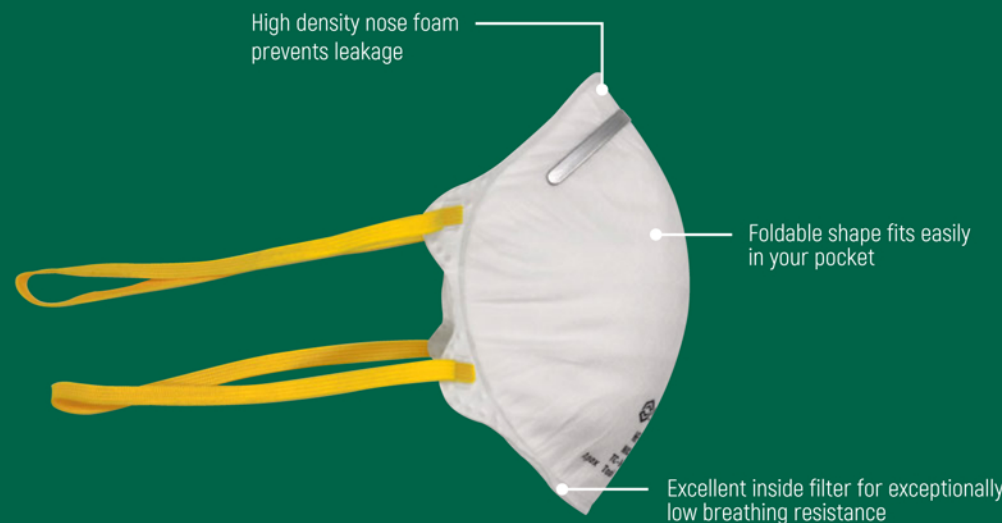

LIMITED WARRANTY
IMPORTANT NOTICE TO PURCHASER:

This limited warranty is made in lieu of the warranties of merchantability, fitness for particular purposes and all other warranties, express or implied. There are no other warranties which extend beyond the description on the face hereof. The physical standards and specifications of Makrite will be met by products sold.

Exclusive Remedies: damages for the breach of this limited warranty are limited to the replacement of such quantity of Makrite products proved to be defectively manufactured. Except as provided above, Makrite shall not be liable or responsible for any loss, damage, or liability, direct, indirect, incidental, special or consequential, arising out of sale, use, misuse or the inability to use products by the user.

WARNINGS:

1. Wearer must read fitting instructions and use limitations before use.

2. This respirator does not protect against the risk of contracting disease or infection.

3. Individuals highly sensitive to natural rubber latex may have an allergic reaction.

4. Before occupational use of this respirator, a written respiratory protection program must be implemented meeting all the local government requirements. In the United States employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing.



RESPIRATOR FITTING INSTRUCTIONS


- 1


- 2


- 3


- 4


- 5



STEP 1

· Hold the respirator in hand with the nosepiece at your fingertips, allowing the headbands to hang freely below your hand

STEP 2

· Press the respirator firmly against your face with the nosepiece on the bridge of your nose

STEP 3

· Stretch and position the top band high on the back of the head. Stretch the bottom band over the head and position below your ears.

STEP 4

· Using both hands, mold nosepiece to the shape of your nose

STEP 5

· To test fit: a) Cup both hands over the respirator being careful not to disturb position, and b) inhale vigorously. If air leaks around the edges, reposition the straps or adjust strap tension for better fit. Please carefully follow these fitting instructions during each use to achieve proper fit.

Executive Standard: USA NIOSH 42 CFR Part 84 - 1995

LOT NO. / MFG Date:

Expiration Date: