

150 MM

20 Disposable Respirators

N95-96K Disposable Respirator

LATEX FREE

Foldable Design



N95-96K



LATEX FREE

Foldable Design


N95 Disposable Respirator

20 Disposable Respirators

LOT NO. / MFG Date:

Expiration Date:

130 MM



Ho Cheng Safety Enterprise Co., Ltd.
No. 5, Alley 9, Lane 48,
San Chun St., Shu Lin District,
New Taipei City 23865, Taiwan
Tel: +886 2 2689 6015

This respirator is approved only in the following configuration:

TC -	Protection ¹	Respirator	Cautions and Limitations ²
84A-7867	N95	N95-96K	ABCJMNOPS

1. Protection

N95 - Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil; time use restrictions may apply.

2. Cautions and Limitations

A - Not for use in atmospheres containing less than 19.5% oxygen.
B - Not for use in atmospheres immediately dangerous to life or health.
C - Do not exceed maximum use concentrations established by regulatory standards.
J - Failure to proper use and maintain this product could result in injury or death.
M - All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA and other applicable regulations.
N - Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
O - Refer to user's instructions, and/or maintenance manuals for information on use and maintenance of these respirators.
P - NIOSH does not evaluate respirators for use as surgical masks.
S - Special or critical user's instructions and/or specific use limitation apply. Refer to user's instruction before donning

¹ - Special or Critical User's Instructions : The model N95-96K filtering facepiece respirator has been manufactured by Approval Holder Makrite for Company Ho Cheng Safety Enterprise Co., Ltd. under TC-84A-7867

SERVICE TIME LIIMITATION

1. Do not use the same respirator for more than 8 hours, continuous or intermittent, in dirty workplaces that could result in high filter loading.
2. Leave contaminated area and remove respirator if breathing becomes difficult or proper fit cannot be obtained or the respirator is damaged, soiled, or distorted. Do not use the same respirator again once it has been used in contaminated area.
3. The maximum service time of this respirator is five (5) consecutive calendar days (including days of non-use), beginning from the first day of use.

130 MM

150 MM



LIMITED WARRANTY

IMPORTANT NOTICE TO PURCHASER:

This limited warranty is made in lieu of the warranties of merchantability, fitness for particular purposes and all other warranties, express or implied. There are no other warranties which extend beyond the description on the face hereof. The physical standards and specifications of Makrite will be met by products sold.

Exclusive Remedies: damages for the breach of this limited warranty are limited to the replacement of such quantity of Makrite, products proved to be defectively manufactured. Except as provided above, Makrite shall not be liable or responsible for any loss, damage, or liability, direct, indirect, incidental, special, or consequential, arising out of sale, use, or misuse, or the inability to use products by user.


Ho Cheng Safety Enterprise Co., Ltd.
No. 5, Alley 9, Lane 48, San Chun St., Shu Lin District, New Taipei City 23865, Taiwan
Tel: 886 2 2689 6015
www.hc-safety.com

Manufactured by Dongguan Sure_On Plastic & Fiber Products Co., Ltd.
ZhuFo Ling Industrial Area, TangXia Town, DongGuan City, GuangDong, China

Made In China

WARNINGS

- Wearer must read fitting instructions and use limitations before use.
- The respirator does not protect against the risk of contracting disease or infection.
- Before occupational use of this respirator, a written respiratory protection program must be implemented meeting all the local government requirements. In the United States employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing.

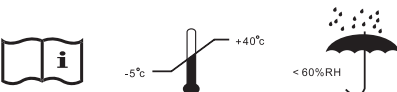

RESPIRATOR FITTING INSTRUCTION



STEP 1

Hold the respirator in hand with the nosepiece at your fingertips, allowing the headbands to hang freely below your hand.



STEP 2

Press the respirator firmly against your face with the nosepiece on the bridge of your nose.



STEP 3

Stretch and position the top band high on the back of the head. Stretch the bottom band over the head and position below your ears.



STEP 4

Using both hands, mold nosepiece to the shape of your nose.



STEP 5

To test fit: a) Cup both hands over the respirator being careful not to disturb position, and b) inhale vigorously. If air leaks around the edges, reposition the straps or adjust strap tension for better fit.

Caution: The respirator should be fit checked prior to each wearing.

175 MM