

220mm

150mm

130mm



## PARTICULATE RESPIRATOR




15  
PCS

ADM102

- Adjustable Nosepiece
- Comfortable Nose Foam
- Easy Breathing & Speaking
- Contour Fit & Foldable Design



### FITTING INSTRUCTIONS

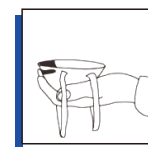


STEP 1 Hold the respirator in hand with the nosepiece at your fingertips, allowing the headbands to hang freely below your hand.




STEP 2 Press the respirator firmly against your face with the nosepiece on the bridge of your nose.



STEP 3 Stretch and position the top band high on the back of the head. Stretch the bottom band over the head and position below your ears.

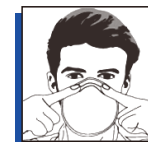


STEP 4 Using both hands' fingertip, mold nosepiece to the shape of your nose.

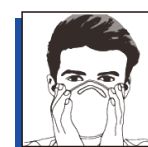

STEP 5 Cup both hands over the particulate respirator and inhale vigorously, a negative pressure should be felt, if not, please re-mold nosepiece following step 4, and/or readjust straps until it fits properly  
Warning: If cannot achieve good seal, do not enter contaminated area.

**NOTE: DO NOT wear the respirator with beards or other facial hair that interferes with direct contact between the face and the edge of the filter.**

ADM102 | 15pcs

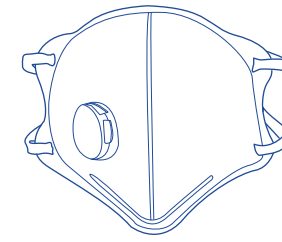

## PARTICULATE RESPIRATOR



New Fortune Co., Ltd  
Room 826, 8F, Three (Southern) Block  
Shushan District Industrial Park Innovation Base, Hefei, Anhui, China.  
TEL: +86-551-62832168 FAX: +86-551-61832199



This respirator is approved only in the following configuration:

TC-	Protection <sup>1</sup>	Respirator	Cautions and Limitations <sup>2</sup>
84A-7869	N95	ADM102	ARC/JMK/OPS

#### 1. Protection

N95-Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil time use restrictions may apply

#### 2. Cautions and Limitations

- A - Not for use in atmospheres containing less than 19.5% oxygen.
- B - Not for use in atmospheres immediately dangerous to life or health.
- C - Do not exceed maximum use concentrations established by regulatory standards.
- J - Failure to properly use and maintain this product could result in injury or death.
- M - All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA and other applicable regulations.
- N - Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
- O - Refer to User's Instructions, and/or maintenance manuals for information on use and maintenance of these respirators.
- P - NIOSH does not evaluate respirators for use as surgical masks.
- S - Special or critical User's Instructions and/or specific use limitations apply. Refer to User's Instructions before donning.

S - Special or Critical User's Instructions: The model ADM102 filtering facepiece respirator has been manufactured by Approval Holder Makrite for New Fortune Co., Ltd. under TC-84A-7869.

ADM102

- Adjustable Nosepiece
- Comfortable Nose Foam
- Easy Breathing & Speaking
- Contour Fit & Foldable Design

15  
PCS



www.goldenscute.com

### USE FOR:

These respirators are suitable with the case of concentration up to 10 times the PEL (Permissible Exposure Limit). They are designed to protect against particles such as those from grinding, sanding, sweeping, sawing, bagging, or coal, iron ore, flour, metal, wood, pollen, and certain other solid substance free of oil. From additional information on the use recommendations for these masks please consult your supervisor or the manufacturer.

1. Failure to follow all instructions and limitations concerning the use of the respirator and/or failure to wear the respirator for the whole duration of exposure to contaminants can seriously reduce the performance characteristics of the respirator and lead to illness, injury or death.

2. Before occupational use of this respirator a written respiratory protection program must be implemented meeting all the local government requirements. In the United States employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing.


3. OSHA standards including 29CFR1910.134 (e) (5) and CSA standard Z94.4-93 require that the wearer be fit-tested.

4. Discard the respirator and replace with a new one if: *Excessive clogging of the respirator causes breathing difficulty. The respirator becomes damaged.*

5. Leave the contaminated area if dizziness, irritation or other distress occurs.

### USE LIMITATIONS:

DO NOT use the respirator or enter or stay in a contaminated area under the following circumstances:  
1. Atmosphere contains less than 19.5% oxygen  
2. For protection against gases or vapors  
3. Contaminants or their concentrations are unknown or immediately dangerous to life or health  
4. Concentrations or contaminants exceed maximum use concentrations in applicable OSHA standards or applicable government regulations or 10 times the PEL (Permissible Exposure Limit), whichever is lower.

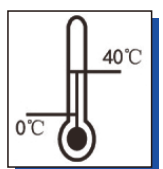

LOT #

Expiration Date: