


USE INSTRUCTIONS

1. Failure to follow all instructions and limitations can seriously reduce the effectiveness of this respirator and could lead to illness, injury or death.
2. A properly selected respirator is essential. Before occupational use, the wearer must be trained by the employer in the correct use of the respirator in accordance with applicable safety and health standards.
3. This product does not supply oxygen. Use only in adequately ventilated areas containing sufficient oxygen to support life.
4. Discard the respirator and replace with a new one if:
 A. excessive clogging of the respirator causes breathing difficulty
 B. the respirator becomes damaged.
5. Leave the contaminated area if dizziness, irritation or other distress occurs.
6. Keep respirators in the display box away from direct sunlight until use.
7. Before occupational use of this respirator, a written respiratory protection program must be implemented meeting all the local government requirements. In the United States, employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing.
8. Before wearing the respirator users must check whether there is crack or spot around the welding spot .if yes,this respirator must be discarded at once and can not be used again.
9. The product is not made with natural rubber latex.
10. S – Special or Critical User Instructions: This respirator has been approved as a NIOSH N95 filtering facepiece respirator, for use in healthcare settings, as a Surgical N95 Respirator conforming to recognized standards for biocompatibility, flammability, and fluid resistance.
Model S-C288S N95 filtering facepiece respirator was manufactured by Approval Holder Guangzhou Harley Commodity Company Limited for Snowcom Medical Commodity (Guangzhou) Co., Ltd under TC-84A-9380.
11. Meet level 3 under 160 mmHg pressure for fluid resistance according to ASTM F 1862.

Donning & Fitting Instructions

- 1.Open the respirator with both hands
- 2.Cup the respirator firmly against your face with the nosepiece on the bridge of your nose.
- 3.Stretch and pull the lower headband over the head and position below your ears.Stretch and pull the top headband on the back of your head above your ears
- 4.Press soft metal or plastic nosepiece to conform snugly around your nose
Using both hands press the nose piece to conform snugly around the nose.
Using only one hand may not ensure a uniform seal.
5. Seal check:
 - To test the fit of a respirator without an exhalation valve, cup both hands over the respirator and exhale sharply.
 - To test the fit of a respirator with an exhalation valve, cup both hands over the respirator and inhale sharply. A negative pressure should be detected inside the respirator.
 - If air flow is felt in the nose area, re-adjust/tighten the nose clip.
 - If air flow is felt around the edges of the respirator, re-position the respirator/head strap to achieve a better fit.
6. Change the respirator immediately if breathing becomes difficult or respirator becomes damaged or distorted.
7. Change the respirator if a proper face seal can not be achieved.
8. Careful observance of these instructions is an important step in safe respirator use.

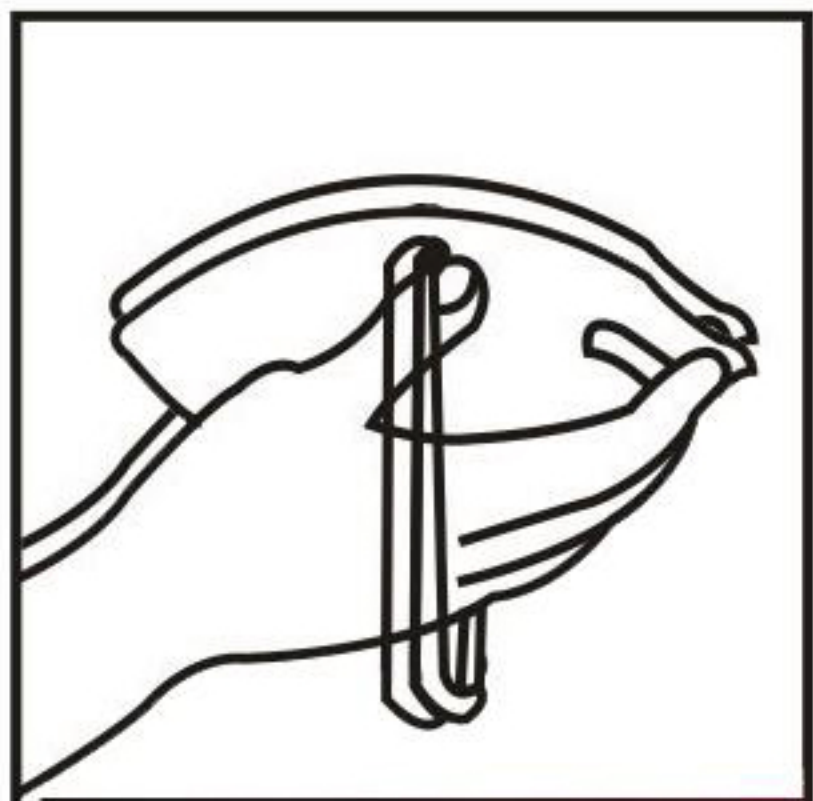

1



2



3



4



Snowcom Medical Commodity (Guangzhou) Co., Ltd

(Plant A01) Yuanzhougang Village,Taiping Town,
Conghua District, Guangzhou, China

TEL: +86-186 2022 5020



THIS RESPIRATOR IS APPROVED ONLY IN THE FOLLOWING CONFIGURATIONS:

TC-	Protection ¹	Respirator	Cautions and Limitations ²
		S-C288S	
84A-9380	Surgical/N95	X	ABCJMNOS

1. Protection:

N95 - Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil.Time restrictions may apply

2. Cautions and Limitations:

- A - Not for use in atmospheres containing less than 19.5 percent oxygen.
- B - Not for use in atmospheres immediately dangerous to life or health.
- C - Do not exceed maximum use concentrations established by regulatory standards.
- J - Failure to properly use and maintain this product could result in injury or death.
- M - All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA, and other applicable regulations.
- N - Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
- O - Refer to User's Instructions, and/or maintenance manuals for information on use and maintenance of these respirators.
- S - Special or Critical User Instruction and/or specific use limitation apply. Refer to User Instructions before donning.