

# SHIELD-95-MEDI

## Surgical N95 Particulate Respirator

User Instructions


IMPORTANT: Keep these user instructions for reference.

248083505-1 R4

### HOW TO WEAR:



1. Turn the respirator so that you are looking at the back side.
2. Pull up on the nose piece and down on the chin pieces to create an opening.
3. Place the respirator on your face so that the center of the nose piece wire is against the bridge of your nose and the chin piece is below your chin.
4. Pull both straps over your head. Move the lower strap around your neck; pull the top strap onto your head.
5. Using one hand, pinch the nose wire across the bridge of your nose and along the sides.
6. Using both hands, press each side of the nose wire strip to conform to your cheek bones on both sides.
7. With the respirator fitted to your face, place finger and thumb on the ends of the Clear View Bar and pinch towards the face to increase downward visibility.

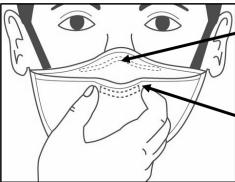

Nose Wire.

The Clearview Bar, once adjusted, should crease in the center and hold the faceplate of the respirator out of user's downward line of sight.

- Seal Check:
- a. Perform a seal check by covering the front of the respirator with both hands, being careful not to disturb the position of the respirator.
  - b. Exhale sharply. A positive pressure inside the respirator means no leakage.
  - c. If leakage occurs, adjust the position and/or tension of the straps and press down along the nose piece as described in previous steps.
  - d. If you cannot achieve a proper seal, remove the respirator and DO NOT enter the contaminated area. See your supervisor.
- For a tighter fit, the bottom head strap can be pulled up around your head.
- To remove respirator; grab the bottom corners with both hands, pull away from your face and upwards so the straps slip off your head.

### ⚠ WARNING

- Before using this respirator, the wearer must read and understand these warnings, instructions, and limitations. Failure to follow them (and any misuse) could reduce the effectiveness of the respirator and result in sickness, serious injury, or death.
1. This respirator helps protect against certain non-oily particulate contaminants but does not eliminate exposure to all particulates or the risk of contracting any disease or infection.
  2. You or your employer, if occupational use, must determine the suitability of this respirator for its intended application. A hazard assessment may determine the need for additional personal protective equipment ("PPE").
  3. Before occupational use of this respirator, a written respiratory protection program must be implemented, meeting all applicable government requirements. In the United States, employers must comply with OSHA 29 CFR 1910.134, which includes medical evaluation, training, and fit testing. In Canada, CSA Standard Z94.4 requirements must be met.
  4. Use only in non-harmful environments.
  5. Keep the respirator in its original packaging before use and always away from contaminants, sunlight, extreme temperatures, and moisture.
  6. Do NOT alter or modify this respirator in any way.
  7. Inspect this respirator before every use. Do NOT use if the respirator is damaged in any way.
  8. Always put on, and wear the respirator in accordance with the user instructions (see package insert).
  9. Make sure that nothing interferes with the seal when you are wearing the respirator. Mustaches, beards or other facial hair, for example, may interfere with the seal and thereby reduce the effectiveness of this respirator.
  10. Evacuate area immediately if your breathing becomes difficult, you experience dizziness, or other respiratory distress occurs.
  11. Do NOT use this respirator in oxygen-deficient environments (less than 19.5% oxygen) or when concentrations of contaminants are immediately dangerous to health or life.
  12. Do NOT reuse respirator. This respirator is for a single use and should be disposed after use.
  13. Do NOT wash the respirator or get it wet. If respirator gets wet, remove and discard
  14. If respirator becomes damaged, soiled, or breathing becomes difficult, dispose of the respirator and use a new one.
  15. Always properly dispose of respirator in accordance with specific local health and safety requirements.
- S – Special or Critical User Instructions: This respirator has been approved as a NIOSH N95 filtering facepiece respirator, for use in healthcare settings, as a Surgical N95 Respirator conforming to recognized standards for biocompatibility, flammability, and fluid resistance at 160 mmHg.

Contact Information:


Lighthouse Worldwide Solutions  
300 West Antelope Road  
White City, OR 97503  
Tel: +1 541-770-5905  
[www.golighthouse.com](http://www.golighthouse.com)

## **MARKING:**



End of shelf-life

YYYY/MM

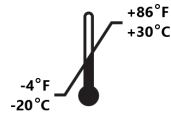

Temperature range



Maximum relative humidity



Lighthouse Worldwide Solutions, Inc.  
White City, OR, USA  
1 (541) 770 5905



**This respirator is approved only in the following configuration:**

TC-	Protection <sup>1</sup>	Respirator	Cautions and Limitations <sup>2</sup>
		SHIELD-95-MEDI	
84A-9431	Surgical/N95	X	ABCJM NOS

### **1. Protection**

N95 – Particulate filter (95% filter efficiency level) effective against particulate aerosols free of oil; time use restrictions may apply

### **2. Cautions and Limitations**

- A Not for use in atmospheres containing less than 19.5 percent oxygen.
- B Not for use in atmospheres immediately dangerous to life or health.
- C Do not exceed maximum use concentrations established by regulatory standards.
- J Failure to properly use and maintain this product could result in injury or death.
- M All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA, and other applicable regulations.
- N Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
- O Refer to User Instructions and/or maintenance manuals for information on use and maintenance of these respirators.
- S Special or critical User's Instructions and/or specific use limitations apply. Refer to User Instructions before donning.