

# GERSON®

## N95 MODEL 3330+ SURGICAL N95 RESPIRATOR

NIOSH Approved TC-84A-9582

Meets CDC Guidelines for TB Exposure Control

This respirator has been approved as a NIOSH N95 filtering facepiece respirator, for use in healthcare settings, as a Surgical N95 Respirator conforming to recognized standards for biocompatibility, flammability, and fluid resistance at 160 mm Hg.

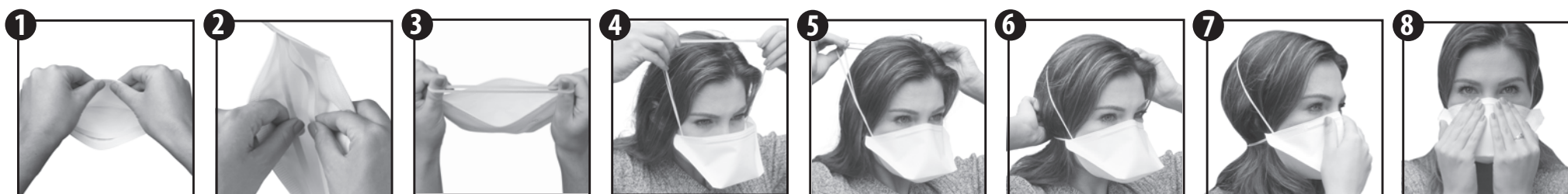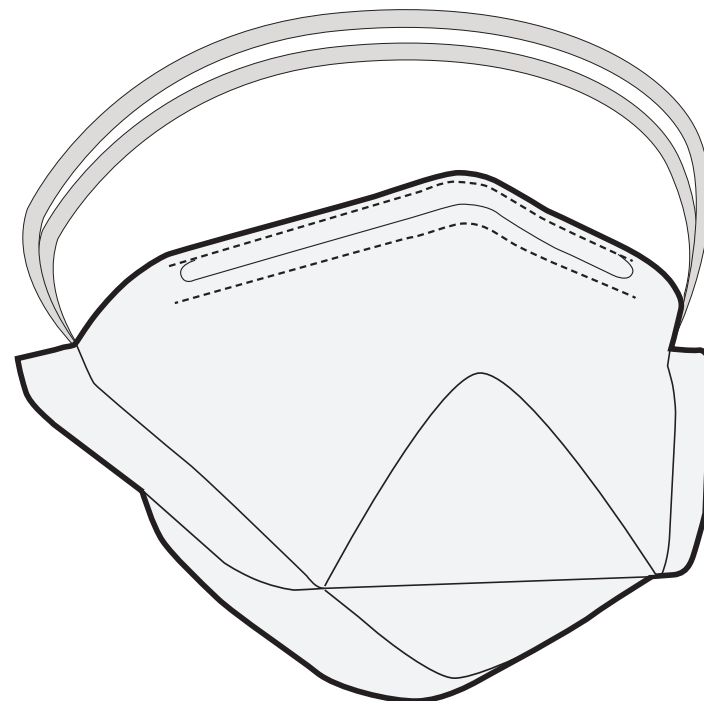1. Unwrap respirator and bend nosepiece into a curve to open respirator.
2. Separate the pre-scored elastic headband into two individual bands.
3. Hold the headbands with your index fingers and thumbs.
4. Cup the open respirator **securely** under your chin.
5. Pull both headbands over your head.
6. Rest the top headband on the crown of your head and bring the bottom one down under your ears.
7. Press the nosepiece around the contour of your nose and cheekbones. Adjust the edges of the respirator to make sure there are no gaps.
8. Perform a Fit Check by gently cupping both hands around the respirator, then exhale sharply to see if you sense or detect any leaks. If leaks are detected, reposition the respirator and/or adjust the straps. If leaks are still detected, then do not enter contaminated area and seek supervisor for further instructions.

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© 2023, All Rights Reserved
Middleboro, MA 02346 U.S.A.
Made in U.S.A.