


# FITTING INSTRUCTIONS

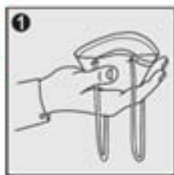

- ① Please the mask on your hand with the nose clip at your fingertips, allowing the head attachment string to hang freely. Then position the lower part of the mask under your chin, and the upper part over the nose.
- ② With the mask over the chin, mouth and the nose, pull the bottom head string over and back of your head and below the ears. Pull the top head string over the top and back your head and above the ears.
- ③ Using both hands, mold the nose area to the shape of your nose by pushing inward while moving your fingertips down both sides of the nose clip.
- ④ To ensure that you have a snug and leak-proof fit, place both hands over the mask and exhale strongly. If the air leaks around the edges, readjust the nose clip and the mask around the edges, and check for air leakage again.