



F333V R95 Particulate Respirator

User's Instructions

WARNING TO USER

Before use, the user must read, understand and follow all User's Instructions, warnings and use limitations of this respirator. This respirator provides limited protection and must be properly selected and fit tested by a competent person to provide protection from known contaminants according to OSHA 29 CFR 1910.134 or applicable safety and health standard. The respirator's face seal must be checked each time before entering a contaminated area and the respirator must be worn during all times of exposure to contaminants. **Failure to follow all of these instructions, warnings and limitations, contained on this package may reduce the effectiveness of this respirator and result in serious sickness or DEATH!!**

The F333V R95 respirator with a filtration efficiency of at least 95% against solid or liquid particles.

DO NOT USE FOR:

The F333V R95 respirator should not be used for gases, vapors, or where concentrations of contaminants exceed OSHA or applicable safety and health standard's maximum use concentrations, are unknown or immediately dangerous to life and health.

USE INSTRUCTIONS:

- Before occupational use of this respirator a written respiratory protection program must be implemented meeting all the local government requirements. In the United States employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing.
- Prior to each use, carefully inspect the entire respirator for correct assembly and for tears or damage including the straps, filter and exhalation valve (if equipped with one). **Discard the entire respirator if any part does not function properly** and get a new respirator.
- If the respirator becomes damaged, distorted, soiled, if breathing becomes difficult or a proper respirator to face seal cannot be maintained, **leave the contaminated area immediately** and get a new respirator.
- If dizziness, irritation, nausea or other distress occurs, **leave the contaminated area immediately** and contact your supervisor or medical professional.
- Do not alter abuse or misuse this respirator.
- Properly dispose of used respirators in accordance with applicable regulations.

USE LIMITATIONS:

- This respirator does not supply oxygen and must not be used in atmospheres containing less than 19.5% oxygen.
- This respirator is not to be used when the user has beards, facial hair, scars or other conditions that prevent a good seal between the face and the respirator.
- This respirator is designed for occupational or professional use by adults who are properly trained in their use and limitations. Respirators must not be used by children.
- Particulates and other contaminants which can be dangerous to your health include those that you cannot see, taste or smell. Respirators can help to protect your lungs against certain airborne contaminants, but cannot prevent entry through other routes such as the skin. Additional types of personal protective equipment may be required.

STORAGE:


- Store unused respirators in their sealed, original packaging in a clean, dry, non-contaminated area that is protected from extreme temperatures, moisture, sunlight, dust and chemicals.
- Never attempt to clean and reuse this respirator.
- Used respirators should not be stored** and should be disposed of properly.



see information supplied by the manufacturer



End of Shelf life 3 Years

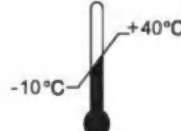

Temperature range of storage conditions -10°C to +40°C



Maximum relative humidity of storage condition <70%

Fido disposable particulate respirator



R95



FILTERING FACEPIECE APPROVAL LABEL

FIDO MASKS CO., LTD.
No. 8-5 She-Ling Li, Yuan-Li Town,
Miao-Li County 35872, Taiwan
TEL: 886-37-741471



THIS RESPIRATOR IS APPROVED ONLY IN THE FOLLOWING CONFIGURATION:

TC-	Protection ¹	Respirator	Cautions and Limitations ²
		F333V	
84A-7795	R95	X	ABCMNOP

1. Protection

R95-Particulate Filter (95% filter efficiency level)
effective against all particulate aerosols;
tune use restrictions may apply

2. Cautions and Limitations

- A - Not for use in atmospheres containing less than 19.5% oxygen.
B - Not for use in atmospheres immediately dangerous to life or health.
C - Do not exceed maximum use concentrations established by regulatory standards.
J - Failure to properly use and maintain this product could result in injury or death.
M - All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA and other applicable regulations.
N - Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
O - Refer to User's Instructions, and/or maintenance manuals for information on use and maintenance of these respirators.
P - NIOSH does not evaluate respirators for use as surgical masks.



FITTING INSTRUCTIONS:

Before fitting the respirator, **carefully inspect the entire respirator** to ensure that the respirator is not damaged and is functioning correctly. **Discard the entire respirator if any part does not function properly.**

- Place the respirator over your nose and mouth. Be sure the nose shape position is on top.
- Pull the bottom strap over your head until it rests below the ears at the back of your neck. Untwist the strap.
- Pull the top strap over your head until it rests on the crown of your head above your ears. Untwist the strap.
- Without removing the respirator, adjust the straps until the respirator is secure. To tighten the respirator, gently pull the ends of both straps. To loosen the respirator, press on the back side of the buckle.
- To check the correct fit, cup with both hands the respirator and exhale powerfully. If any airflow indicates leakage around the edges of the respirator stretch the headbands again by pulling them backwards. Then correct the position of the respirator (see picture) and check again. **WARNING!** It is not allowed to use the mask without proper fit. **NEVER USE THE MASK AND DO NOT ENTER A CONTAMINATED AREA WITHOUT A PROPER SEAL.**

Note: Pinching the nosepiece using one hand may result in improper fit and less effective respirator performance. Use two hands.

Manufacturer Address:

Fido Mask Co., Ltd
8-5 She-Ling, Yuan-Li Town, Miao-Li County 35872, Taiwan

Produce date: YYYY/MM

Lot No. XXXXXX

(the Lot No. XXXXXX will be as same as abbreviated Label indicated, the "XXXXXX" should be printed as purchase order number)

Contents: 10 pcs.

Fido disposable particulate respirator



R95