



Jiangyin Chang-hung Industrial Manufacturing Factory LLC

----Professional PPE Manufacturer----

Address: No.16, Huannan Rd., Huang Tang, Xiake Town, Jiangyin City, Jiang Su, China

TEL: +86-510-86913080

FAX: +86-510-86913082

E-Mail: marketing@grandemask.com

Web: <http://www.grandemask.com>


N95 Single-use Particulate Respirator Donning Process (With Adjustable Buckle)

Check prior to use:

1. Check that there are no cracks or tears around the welds on the mask. If there are cracks or tears the respirator must be disposed and a new mask used.
2. Check the straps are firmed enough to pull over your head.
3. Check the nosepiece and inner sponge are glued closely on the respirator.
4. Check the valve is set up well if the respirator is with exhalation valve.
5. Check the valid date which is printed on the bottom of box.

Fitting instructions:

1. Before use, the wearer must first be trained by the employer in the correct use of the respirator in accordance with applicable safety and health standards for the contaminant and exposure level in the assigned work area. A written respiratory protection program must be implemented, meeting all the requirements of OSHA 29 CFR 1910.134, including training, medical evaluation and fit testing. For more information, please see OSHA's website at <http://www.osha.gov>.
2. Cup the respirator in your hand with the nosepiece at your fingertips allowing the headbands to hang freely below your hand.
3. Cup the respirator firmly against your face with the nosepiece on the bridge of your nose.
4. Pull the top strap over your head resting it high at the top back of your head. Pull the bottom strap over your head and position it around the neck below the ears. Pull strap at either side of buckle to adjust length of strap for various facial sizes.
5. Press soft nosepiece to conform snugly around the nose.
6. CHECK FOR LEAK-TIGHTNESS (FACIAL FIT) AS FOLLOWS:
 - *place both hands over the respirator and exhale sharply
 - *positive pressure should be inside the respirator
 - *if you feel air escaping around the edges, re-adjust the respirator by tightening the nose-piece
 - *repeat the procedure until satisfactory face-fit has been achieved.
7. Careful observance of these instructions is an important step in safe respirator use.

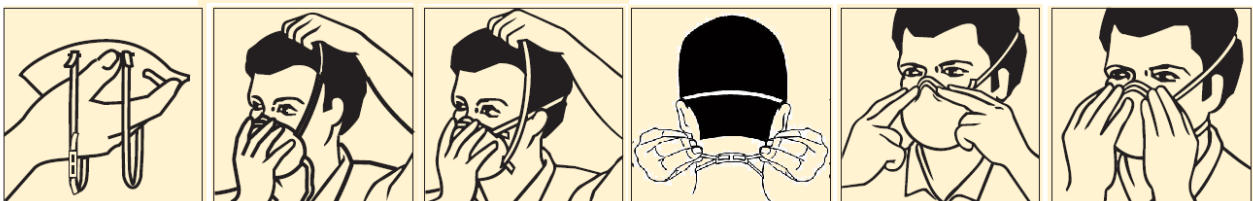



Jiangyin Chang-hung Industrial Manufacturing Factory LLC

----Professional PPE Manufacturer----

Address: No.16, Huannan Rd., Huang Tang, Xiake Town, Jiangyin City, Jiang Su, China

TEL: +86-510-86913080

FAX: +86-510-86913082

E-Mail: marketing@grandemask.com

Web: <http://www.grandemask.com>

Use against:

1. Particulate aerosols free of oil.

Do not use against:

1. Concentrations of contaminants which are unknown or are immediately dangerous to life or health.
2. Concentrations of particulates which exceed the maximum use concentration or 10 times the OSHA Permissible Exposure Limit, whichever is lower.
3. Gases, vapors, asbestos, paint spray, sandblasting, or particulate materials which generate harmful vapors.
4. Oil-based mists.

Warnings:

1. Failure to follow all instructions and limitations on the use of this product could reduce the effectiveness of respirator and result in sickness or death.
2. A properly selected respirator is essential to protect your health. Before using this respirator, consult an industrial hygienist or Occupational Safety Professional to determine the suitability for your intended use.
3. This respirator does not supply oxygen and must not be used in atmospheres containing less than 19.5% oxygen.
4. Do not use when concentrations of contaminants are immediately dangerous to health or life.
5. Leave work area immediately and return to fresh air if breathing becomes difficult, or dizziness or other distress occurs.
6. Do not use with beards or other facial hair that may inhibit contact between the face and the edge of the respirator.
7. Never alter, modify or abuse this respirator.
8. Only for single use. No maintenance necessary. Discard the used respirator after single use.
9. Change respirator immediately if breathing becomes difficult or respirator becomes damaged or distorted, or a proper face fit can't be maintained.
10. Store unused respirators in box in a clean, dry, non-contaminated area.