




**JACKSON\*** 64260  
SAFETY Brand

## **R10** Particulate Respirator with Comfort Strap and Dual Valves N95

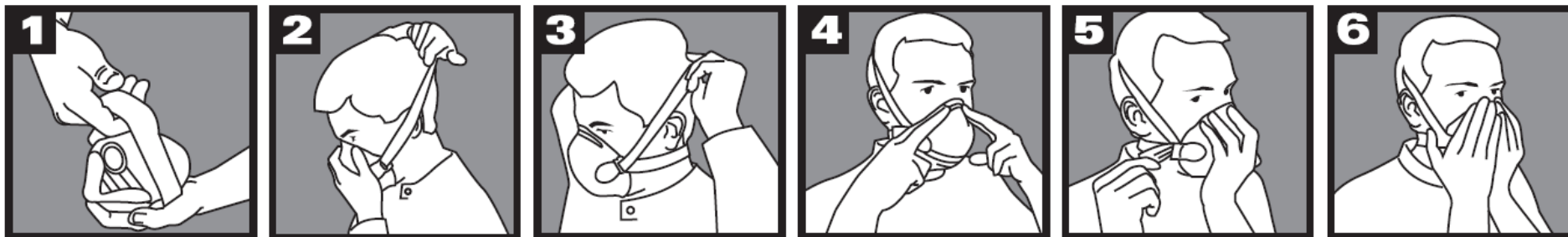

Ⓜ Follow these instructions each time the respirator is worn. (See pictures 1-6)

1. Cup the body of the mask in one hand and grasp both head straps with the other hand. 2. Place the mask over the mouth and nose and pull the head straps over the head. 3. Pull one of the head straps up over the crown of the head and the other behind the neck, taking care to keep the straps flat. 4. Using both hands, carefully mould the nose wire to the shape of your nose. 5. The overall tension can be increased by pulling the tab on the right side buckle. Adjust the tension between the top and bottom headbands by sliding the strap through the left buckle. 6. Test the seal of the mask to the face. Place both hands over the respirator. Inhale and exhale sharply. If any air leaks around the nose, re-adjust the nose wire. If any air leaks around the edges, tighten the head straps and/or reposition the respirator. If you cannot achieve a proper seal, do not enter the contaminated area and see your supervisor. Please notice that you will feel air flow across the lateral valves, do not confuse that with air leaks.