



Kimberly-Clark  
1400 Holcomb Bridge Road  
Roswell, GA 30076-2199  
Tel: (800)331-4112



This respirator is only approved in the following configuration :

TC -	Protection <sup>1</sup>	Respirator		Cautions and Limitations <sup>2</sup>
		NPC001	NPC001V	
84A - XXXX	N95	X		ABCJMNOPS
84A - XXXX	N95		X	ABCJMNOPS

#### 1. Protection

N95 - Particulate Filter ( 95% filter efficiency level ) effective against particulate aerosols free of oil ;  
time use restrictions may apply

#### 2. Cautions and Limitations

- A - Not for use in atmospheres containing less than 19.5% oxygen.
- B - Not for use atmospheres immediately dangerous to life or health.
- C - Do Not exceed maximum use concentrations established by regulatory standards.
- J - Failure to properly use and maintain this product could result in injury or death.
- M - All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA and other applicable regulations.
- N - Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
- O - Refer to users instructions, and / or maintenance manuals for information on use and maintenance of these respirators.
- P - NIOSH does not evaluate respirators for use as surgical masks.
- S - Special or critical User's Instructions and/or specific use limitations apply.  
Refer to User's Instructions before donning.

#### CAUTION :


Filtering face pieces are to be inspected to each use assure there are no holes in the breathing zone other than the punctures around staples and no damage has occurred, Enlarged holes resulting from ripped or torn filter material around staples are considered damage.

#### LIMITED WARRANTY

##### IMPORTANT NOTICE TO PURCHASER:

This limited warranty is made in lieu of the warranties of merchantability, fitness for particular purposes and all other warranties, express or implied. There are no other warranties which extend beyond the description on the face hereof. The physical standards and specifications on NIPPON PAPER CRECIA LTD. will be met by products sold.


Exclusive Remedies: damages for the breach of this limited warranty are limited to the replacement of such quantity of NIPPON PAPER CRECIA LTD., products proved to be defectively manufactured. Except as provided above, NIPPON PAPER CRECIA LTD. shall not be liable or responsible for any loss, damage, or liability, direct, indirect, incidental, special, or consequential, arising out of sale, use, or misuse, or the inability to use products by the user.

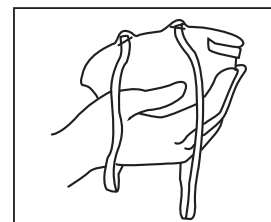

#### FITTING INSTRUCTIONS :

##### Step 1

Hold respirator in your hand as shown, with the nosepiece at your fingertips and the headbands hanging freely below your hand.

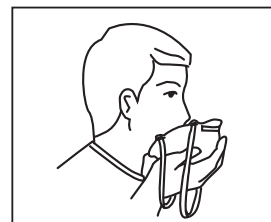

##### Step 2

Place respirator under your chin with the nosepiece up. Pull the top headband to top back of head.



##### Step 3

Pull the bottom headband overhead and place around neck below ears.



##### Step 4

Adjust headbands for comfortable fit. Push the nosepiece inward with your fingertips to conform to the shape of your nose.

#### WARNINGS :

1. Wearer must read filtering instructions and use limitations and inspect the straps before use. If there are any cracks around or in the welded area, then the mask is considered damaged and should be discarded.
2. This respirator does not protect against the risk of contracting disease or infection.
3. Before use, must be fit tested either quantitatively or qualitatively by a qualified fit tester before using the respirator in a contaminated environment, as per the OSHA guidelines, 1910.134.

S - Special or critical User's Instructions: The model NPC001 & NPC001V filtering face piece respirator has been manufactured by Approval Holder Makrite Industries Inc. for Kimberly-Clark Corporation.