N95-NIOSH Approved Particulate Respirator *Nuisance Level OV-AG Relief**

GERSON® 1735

N95-NIOSH Approved Particulate Respirator *Nuisance Level OV-AG Relief**

GERSON®
Leading the Way

20
CUPPED/NESTED

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NOTE: Before issuing any respirator to be worn, a face fit test must be performed. See OSHA regulation 1910.134.

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1. Hold the respirator in hand with nosepiece at your fingertips, allowing the headbands to hang freely below your hand.




2. Press the respirator firmly against your face with the nose-piece on the bridge of your nose.

3. Stretch and position the top band high on the back of your head. Stretch the bottom band over the head and position it below your ears.

4. Using both hands, mold the metal nosepiece to the shape of your nose.

5. To test fit cup both hands over the respirator and exhale vigorously. If air flows around your nose, tighten the nosepiece, if air leaks around the edges, reposition the straps for a better fit.

6. Remember, careful observance of these fitting instructions is an important step in safe respirator use.

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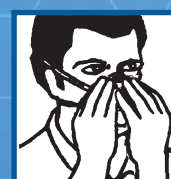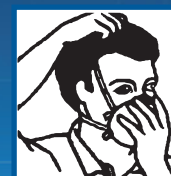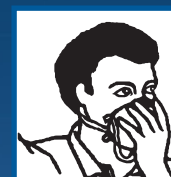* Gerson recommended for relief against the irritating effects of nuisance levels of organic vapor or acid gas. This product does not provide lung protection from organic vapors or acid gases. Nuisance level refers to concentrations not exceeding OSHA PEL or applicable government occupational exposure limits, whichever is lower.