


FITTING INSTRUCTIONS: Thread top strap through buckles. Repeat for bottom strap and buckle.



- 1. Position bottom strap around the back of the head, just below the ears.**

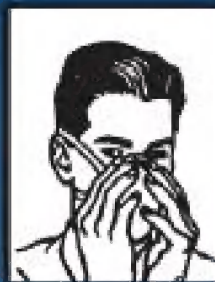

- 4. Using both hands, mold the metal nosepiece to the shape of your nose.**



- 2. Position top strap around the back of the head just above the ears.**



- 5. To release strap tension without removing the respirator, push out on back of buckle.**



- 3. Adjust the tension by pulling on end tab of each strap.**



- 6. Check seal of respirator on the face, prior to use; for proper fit. A) Cup both hands over the front of the respirator being careful not to disturb position. B) Inhale sharply. A negative pressure should be felt inside respirator. If leakage is suspected, adjust position of respirator and/or strap and repeat procedure.**

NOTE: If a proper fit cannot be achieved, do not enter contaminated area, See your supervisor.