


N95 Particulate Respirator Mask

(95% Filtration Efficiency)

MAGNUM




Single Use / Latex Free

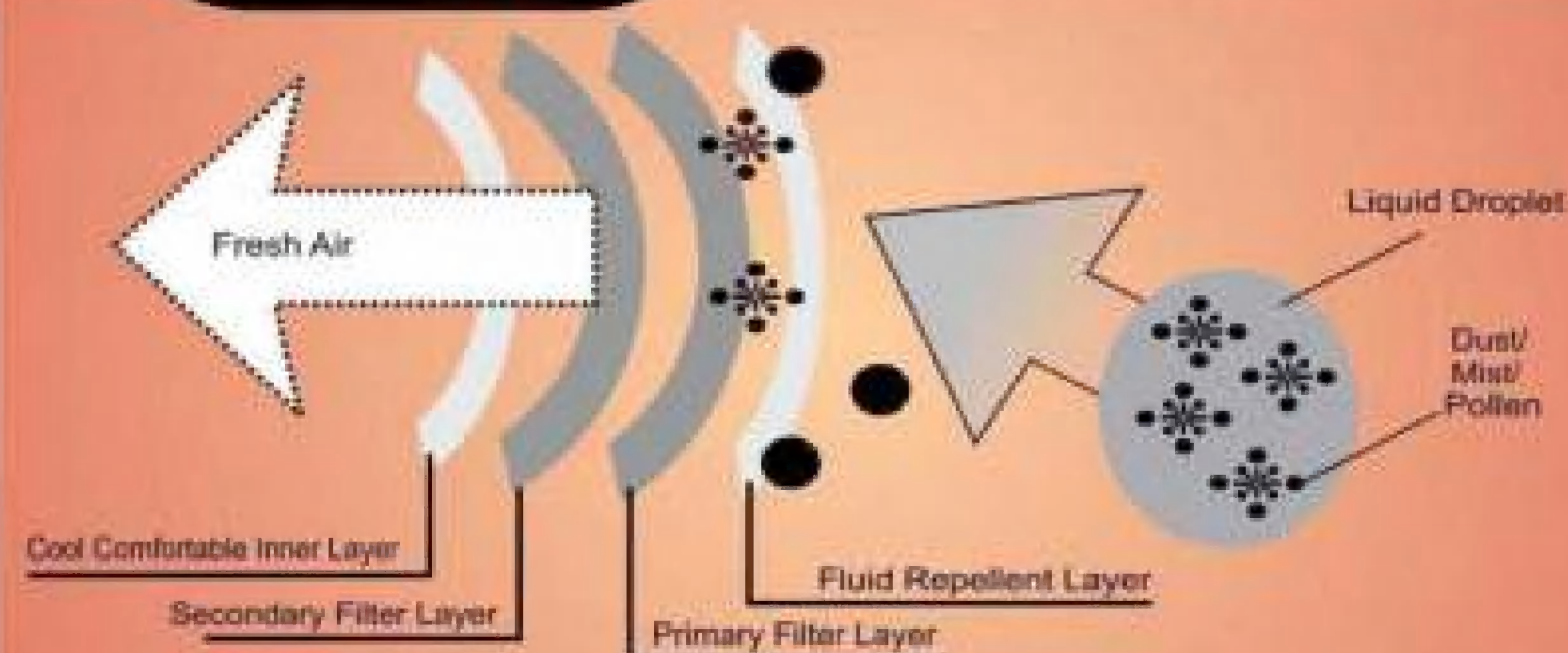

User Instructions :



1) Spot the two head straps inside mask and pull it away.

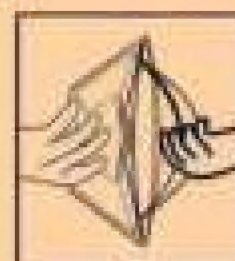

2) Open the mask and ensure nose clip facing upward.



3) Place the respirator under your chin and pull the top headband over your head and above the ears. Pull the bottom elastic and place it on your neck.



4) Press adjustable nose clip to prevent leakage.



5) Push the nose clip inward with your finger tip to confirm to the shape of your nose.



6) To check fit, place both hands completely over the filter mask and exhale. If air leaks around nose adjust the nose clip.

Manufactured by :

MAGNUM MEDICARE PVT. LTD.

Survey No. 14, Hissa No.2, Unit No. 21,22,23 Supreme Estate, Village Devdal, Taluka Vasai, Thane - 401203. Tel. : 022-28484282 • Customer Care No. : 08080850482
e-mail : info@magnummedicare.com

MAGNUM

Warnings :

- 1) A properly selected respirator is essential to protect your health. Before using this or any product, consult an Industrial Hygienist or Occupational Professional to determine its suitability for your intended use.
- 2) Reference OSHA respiratory protection standard 29 CFR 1910.134.
- 3) This product is for dry aerosols only. It will not protect against gases or vapours or in spray paint operations.
- 4) Do not use in sand blasting operations.
- 5) Do not use for protection against asbestos.
- 6) This product does not supply oxygen. Use only in well ventilated areas containing sufficient oxygen to support life.
- 7) Do not use when concentrations of contaminants are unknown or immediately dangerous to life, or when concentrations exceed 10X PEL.
- 8) Leave work area immediately if: (a) breathing becomes difficult; (b) dizziness or other signs of distress occur.
- 9) If user has any illness or disorders, consult a licensed medical doctor to determine suitability of use.
- 10) Facial hair and certain facial characteristics may prohibit effective use of the product.
- 11) Use only in accordance with instructions relating to this product.
- 12) Discard and replace respirator if it becomes damaged or breathing resistance becomes excessive.
- 13) Do not alter or modify this product in any way.
- 14) Follow the manufacturer's instructions for fitting the respirator to the face.

NOTE : Before issuing any respirator to be worn, a face fit test must be performed.



MAGNUM MEDICARE PVT LTD



Correspondence: 602 Manta Heights, Devidas Lane, Borivali(West) Mumbai-400 103, Maharashtra, India

Factory: Survey No.14, Hissa No.2, Unit No.21,22,23 Supreme Estate,Village Devdal, Taluka Vasai,Thane

Tel:022-65948415. Tel/fax:022-28484282 E-mail:info@magnummedicare.com or magnum.dsp@gmail.com

These respirators are approved only in the following configurations:

TC-	Protection ¹	Respirator	Cautions and limitations ²
		N95001	
84A-N95001a	N95	X	ABCDMNOP

1. Protection

N95-Particulate Filter (95% filter efficiency level) effective against particulate - aerosols
free of oil; time use restrictions may apply

2.Cautions and Limitations

A-Not for use in atmospheres containing less than 19.5 percent oxygen.

B-Not for use in atmospheres immediately dangerous to life or health.

C-Do not exceed maximum use concentrations established by regulatory standards.

J-Failure to properly use and maintain this product could result in injury or death.

M-All approved respirators shall be selected, fitted, used and maintained in accordance with MSHA, OSHA, and other applicable regulations.

N-Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.

O-Refer to User's Instructions and / or maintenance manuals for information on use and maintenance of these respirators.

P-NIOSH does not evaluate respirators for use as surgical masks.

Manufactured by :

MAGNUM MEDICARE PVT. LTD.

Survey No. 14, Hissa No.2, Unit No. 21,22,23 Supreme Estate, Village Devdal, Taluka Vasai, Thane - 401203. Tel. : 022-28484282 • Customer Care No. : 08080850482
e-mail : info@magnummedicare.com