


FITTING INSTRUCTIONS:


Step 1
Hold respirator in your hand as shown, with the nosepiece at your fingertips and the headbands hanging freely below your hand.



Step 2

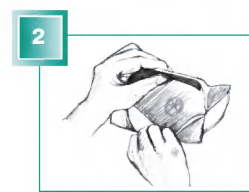

Step 3
To fit check: a) cup both hands over the mask being worn; b) inhale and exhale vigorously. If air leaks around the edges during inhalation or exhalation, re-position the straps or adjust strap tension to better fit. Repeat (a), (b) and (c). If there are still noticeable leaks then this mask may not be suitable for you and we advise you to seek another style from our range.

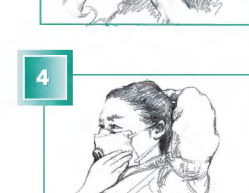

Step 4
Careful observance of these fitting instructions is an important step in safe respirator use.

Caution:
The respirator should be fit checked for each wearing. Place both hands completely over the respirator and exhale. If air leaks around nose or the edges, adjust the nosepiece and/or headbands until a good fit is achieved.