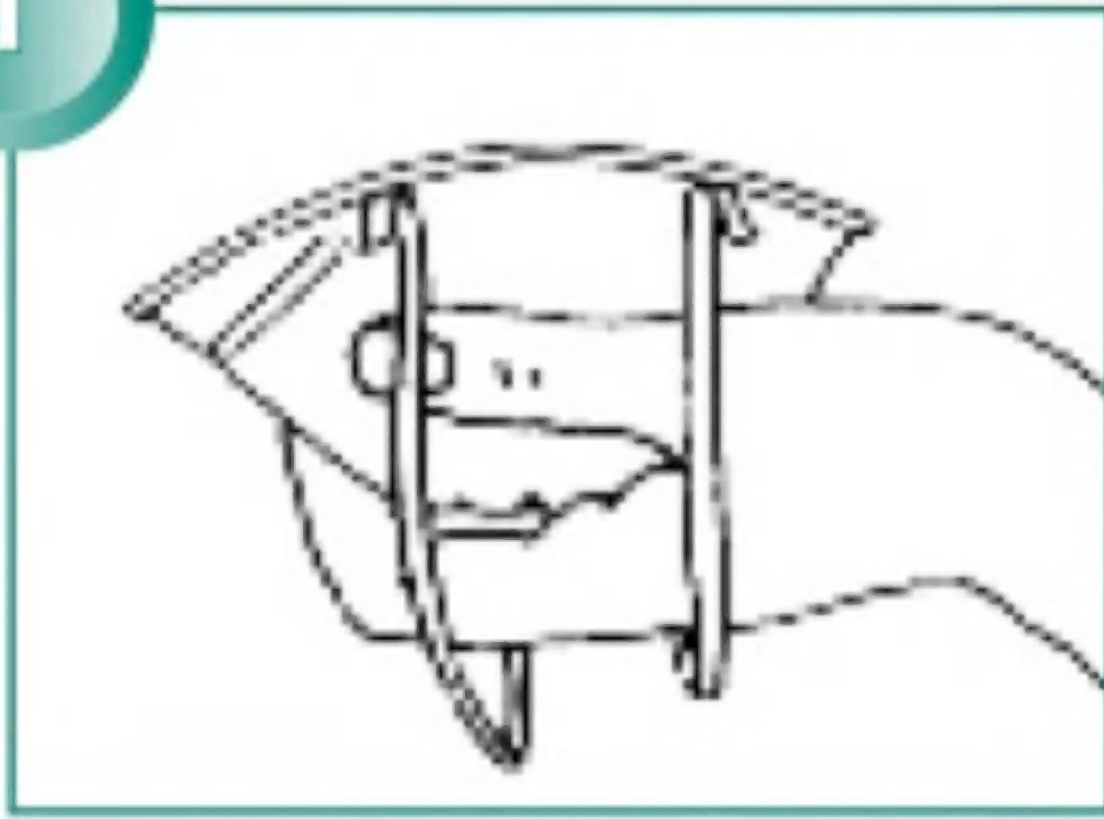


**1**

## **FITTING INSTRUCTIONS:**

### **Step 1**

Hold respirator in your hand as shown, with the nosepiece at your fingertips and the headbands hanging freely below your hand.

**2**

### **Step 2**

Place respirator under your chin with the nosepiece up. Pull the top headband to top back of head.

**3**

### **Step 3**

Pull the bottom headband overhead and place around neck below ears.

**4**

### **Step 4**

Adjust headbands for comfortable fit. Push the nosepiece inward with your fingertips to conform to the shape of your nose.

### **Caution:**

The respirator should be fit checked for each wearing. Place both hands completely over the respirator and exhale. If air leaks around nose or the edges, adjust the nosepiece and/or headbands until a good fit is achieved.