

140 MM

20  
Particulate  
Respirators

N95  
Particulate  
Respirators

N95-001

Innovative Healthcare Solutions  
MedaCure®

MedaCure®  
Innovative Healthcare Solutions

N95-001

NIOSH  
Approved  
Product



N95  
Particulate  
Respirators


20  
Particulate  
Respirators

130 MM


140 MM

130 MM

200 MM



**MEDACURE INC.**  
545 N. Arlington Ave.,  
East Orange, NJ 07017 USA  
T: 718-596-1120  
www.medacureinc.com



This respirator is approved only in the following configuration:

TC-	Protection <sup>1</sup>	Respirator	Cautions and Limitations <sup>2</sup>
84A-5411	N95	N95-001	X ABCJMNOPS

1. Protection

N95-Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil-time use restrictions may apply

2. Cautions and Limitations

A - Not for use in atmospheres containing less than 19.5% oxygen.  
B - Not for use in atmospheres immediately dangerous to life or health.  
C - Do not exceed maximum use concentrations established by regulatory standards  
J - Failure to properly use and maintain this product could result in injury or death.  
M - All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA and other applicable regulations  
N - Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.  
O - Refer to user's instructions, and/or maintenance manuals for information on use and maintenance of these respirators.  
P - NIOSH does not evaluate respirators for use as surgical masks.  
S - Special or critical User's Instructions and/or specific use limitations apply. Refer to User's Instructions before donning.

S - Special or Critical User's Instructions: The model N95-001 filtering face piece respirator has been manufactured by Approval Holder Makrite for Medacure under TC-84A-5411.

**SERVICE TIME LIMITATION**

- Do not use the same respirator for more than 8 hours, continuous or intermittent, in dirty workplaces that could result in high filter loading.
- Leave contaminated area and remove respirator if breathing becomes difficult or proper fit cannot be obtained or the respirator is damaged, soiled, or distorted. Do not use the same respirator again once removed.
- The maximum service time of this respirator is five (5) consecutive calendar days (including days of non-use), beginning from the first day of use.

**MEDACURE INC.**  
545 N. Arlington Ave. East Orange, NJ 07017 USA  
T: 718-596-1120  
www.medacureinc.com

MedaCure®  
Innovative Healthcare Solutions

N95-001

NIOSH  
Approved  
Product



N95  
Particulate  
Respirators

20  
Particulate  
Respirators

**LIMITED WARRANTY  
IMPORTANT NOTICE TO  
PURCHASER:**


This limited warranty is made in lieu of the warranties of merchantability, fitness for particular purposes and all other warranties, express or implied. There are no other warranties which extend beyond the description on the face hereof. The physical standards and specifications of Medacure will be met by products sold.

Exclusive Remedies: damages for the breach of this limited warranty are limited to the replacement of such quantity of Medacure products proved to be defectively manufactured. Except as provided above, Medacure shall not be liable or responsible for any loss, damage, or liability, direct, indirect, incidental, special or consequential, arising out of sale, use, misuse or the inability to use products by the user.

Made in China

**WARNINGS**

- Wearer must read fitting instructions and use limitations before use.
- This respirator does not protect against the risk of contracting disease or infection.
- Before occupational use of this respirator, a written respiratory protection program must be implemented meeting all the local government requirements. In the United States employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing.

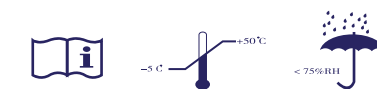

**RESPIRATOR FITTING INSTRUCTION:**



**STEP 1**  
Hold the respirator in hand with the nosepiece at your fingertips, allowing the headbands to hang freely below your hand.



**STEP 2**  
Press the respirator firmly against your face with the nosepiece on the bridge of your nose.



**STEP 3**  
Stretch and position the top band high on the back of the head. Stretch the bottom band over the head and position below your ears.



**STEP 4**  
Using both hands, mold nosepiece to the shape of your nose.



**STEP 5**  
To test fit: a) Cup both hands over the respirator being careful not to disturb position, and b) inhale vigorously. If air leaks around the edges, reposition the straps or adjust strap tension for better fit. Please carefully follow these fitting instructions during each use to achieve proper fit.

LOT #

Expiration Date: