MAKRITE INDUSTRIES INC.  
11F-5, No. 79 Sec. 1, Hsin Tai Wu Rd.  
Hsichih, Taipei Hsien, Taiwan  
Tel: (02)26982419

This respirator is approved only in the following configuration:

TC -	Protection <sup>1</sup>	Respirator	Cautions and Limitations <sup>2</sup>
84A- 6253	N95	MK910-N95DS2 x	ABCJMNOP

1. Protection

N95-Particulate Filter (95% filter efficiency level)effective against particulate aerosols free of oil;time use restrictions may apply.

2. Cautions and Limitations

- A - Not for use in atmospheres containing less than 19.5% oxygen.
- B - Not for use in atmospheres immediately dangerous to life or health.
- C - Do Not exceed maximum use concentrations established by regulatory standards
- J - Failure to properly use and maintain this product could result in injury or death.
- M - All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA and other applicable regulations
- N - Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
- O - Refer to users instructions, and/or maintenance manuals for information on use and maintenance of these respirators.
- P - NIOSH does not evaluate respirators for use as surgical masks.

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
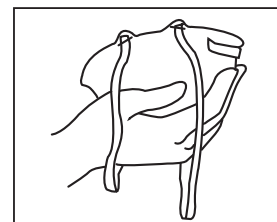
**FITTING INSTRUCTIONS :**

**Step 1**

Hold respirator in your hand as shown, with the nosepiece at your fingertips and the headbands hanging freely below your hand.

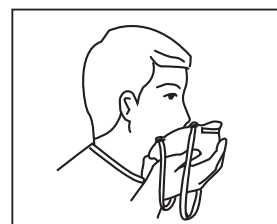**Step 2**

Place respirator under your chin with the nosepiece up. Pull the top headband to top back of head.


**Step 3**

Pull the bottom headband overhead and place around neck below ears.

**Step 4**

Adjust headbands for comfortable fit. Push the nosepiece inward with your fingertips to conform to the shape of your nose.

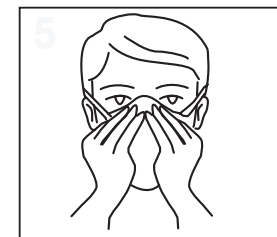**Step 5**

Cup both hands over the mask, checking for tightness, exhale sharply. If air leaks around the edges, adjust the nosepiece and/or headbands until a good fit is achieved.

**WARNINGS:**

1. Wearer must read fitting instructions and use limitations before use.
2. This respirator does not protect against the risk of contracting disease or infection.
3. Individuals highly sensitive to natural rubber latex may have an allergic reaction.
4. Before use, the wearer must be trained in the proper use, including fit testing, in accordance with applicable safety and health standards for the containment and exposure level in the work area as per the OSHA guidelines, 29 CFR 1910.134.

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