



MAKRITE INDUSTRIES INC.
11F-5, No. 79 Sec. 1, Hsin Tai Wu Rd.
Hsichih, Taipei Hsien, Taiwan
Tel: (02)26982419



This respirator is approved only in the following configuration:

TC –	Protection ¹	Respirator	Cautions and Limitations ²
84A- 6270	N95	MK910-N95DS2V x	ABCJMNOP

1. Protection

N95-Particulate Filter (95% filter efficiency level)effective against particulate aerosols free of oil;time use restrictions may apply .

2. Cautions and Limitations

- A - Not for use in atmospheres containing less than 19.5% oxygen.
- B - Not for use in atmospheres immediately dangerous to life or health.
- C - Do Not exceed maximum use concentrations established by regulatory standards.
- J - Failure to properly use and maintain this product could result in injury or death.
- M - All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA and other applicable regulations.
- N - Never substitute, modify, add, or omit parts. Use only exect replacement parts in the configuration as specified by the manufacturer.
- O - Refer to User's Instructions, and/or mainterance manuals for information on use and maintenance of these respirators.
- P - NIOSH does not evaluate respirators for use as surgical masks.

このマスクは次の形式にのみ認められています


T C –	保護 ¹	マスク	警告と制限 ²
		MK910-N95DS2V	
84A - 6270	N95	x	ABCJMNOP

1.保 護

N95 – 油分を含まない粒子エアロゾルに対して有効な粒子フィルター(95%フィルター効果レベル) 使用時間制限が適用される可能性があります。

2. 注意事項および制限


- A- 酸素濃度 19.5%未満の環境では使用できません。
- B- 生命または健康に対して直接的な危険性のある環境では使用できません。
- C- 規格基準により設定された最大使用濃度を超えてはいけません。
- J- 本製品を正しく使用または維持しなかった場合は怪我または死に至る可能性があります。
- M- 全ての認可された防じんマスクはMSHA, OSHA および他の適用規格に従い、選択・着用・使用 され、そして維持されなければなりません。
- N- 決して部品を代用したり、改造・追加・省略などしてはなりません。
メーカーにより指定された構成の中の正しい交換部品のみを使用してください。
- O- これらの防じんマスクの使用とメンテナンスに関する情報は、取扱説明書もしくはメンテナンスマニュアルを参照してください。
- P- NIOSH は防じんマスクを手術用マスクとして使用することは勧めません。

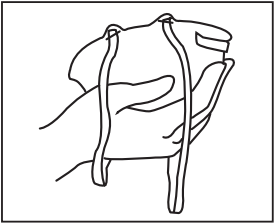

FITTING INSTRUCTIONS :


Step 1

Hold respirator in your hand as shown, with the nosepiece at your fingertips and the headbands hanging freely below your hand.

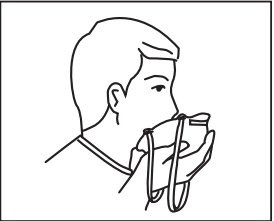

Step 2

Place respirator under your chin with the nosepiece up. Pull the top headband to top back of head.

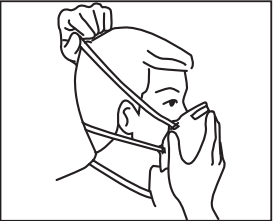

Step 3

Pull the bottom headband overhead and place around neck below ears.



Step 4

Adjust headbands for comfortable fit. Push the nosepiece inward with your fingertips to conform to the shape of your nose.



Step 5

Cup both hands over the mask, checking for tightness, exhale sharply. If air leaks around the edges, adjust the nosepiece and/or headbands until a good fit is achieved.

WARNINGS:

1. Wearer must read fitting instructions and use limitations before use.
2. This respirator does not protect against the risk of contracting disease or infection.
3. Individuals highly sensitive to natural rubber latex may have an allergic reaction.
4. Before use, the wearer must be trained in the proper use, including fit testing, in accordance with applicable safety and health standards for the containment and exposure level in the work area as per the OSHA guidelines, 29 CFR 1910.134.