


FITTING INSTRUCTIONS:

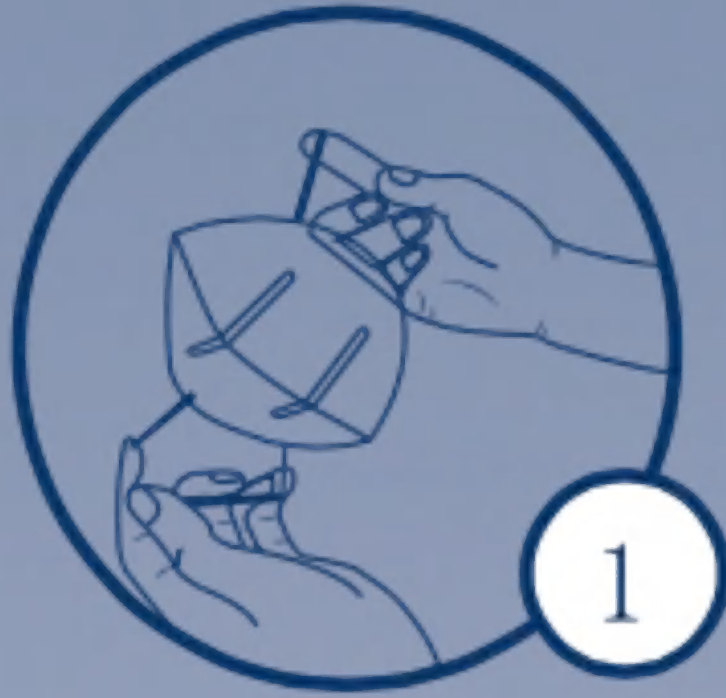

Hold the respirator in your hand as shown, with the left hand holding the left strap and the right holding the right strap.

Place the respirator over the nose, mouth and chin. Pull the left and right strap so that the top of the strap goes around the ear. Hold the hook attached to the right strap and hook onto the left strap on the back of the head for a proper

To test it: Cup both hands over the mask. Inhale sharply. Negative pressure should be felt inside the mask. If air leaks around the edges, reposition the straps or adjust the mask and hook for a better fit.

Caution:

The respirator should be fit checked prior to each wearing.