


FITTING INSTRUCTIONS



STEP 1 Hold the respirator in hand with the nosepiece at your fingertips, allowing the headbands to hang freely below your hand.

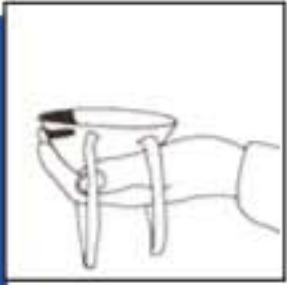

STEP 2 Press the respirator firmly against your face with the nosepiece on the bridge of your nose.



STEP 3 Stretch and position the top band high on the back of the head. Stretch the bottom band over the head and position below your ears.



STEP 4 Using both hands' fingertip, mold nosepiece to the shape of your nose.



STEP 5 Cup both hands over the particulate respirator and inhale vigorously, a negative pressure should be felt, if not, please re-mold nosepiece following step 4, and/or readjust straps until it fits properly

Warning: If cannot achieve good seal, do not enter contaminated area.

NOTE: DO NOT wear the respirator with beards or other facial hair that interferes with direct contact between the face and the edge of the filter.