

140 MM

130 MM


130 MM

140 MM


Silicone Rubber Seal Respirator N95  
95% Efficient Protection Against  
Non-Oil Particulates

9800-N95  
TC-84A-XXXX  
20pcs/pack





Makrite Industries Inc.  
11F-S, No.79, Sec. 1,  
Xintai 5th Rd., Xizhi Dist.,  
New Taipei City 221, Taiwan  
Tel: +886-2-2698-2419



This respirator is approved only in the following configuration:

TC- Protection	Respirator 9800-N95	Silicone Seal S9800	Cautions and Limitations <sup>2</sup>
84A-XXXX	N95 x	x	ABCJMNOPS

1. Protection

N95-Particulate Filter (95%  
filter efficiency level)effective  
against particulate aerosols  
free of oil; use restrictions  
may apply

2. Cautions and Limitations

A - Not for use in atmospheres containing less than 19.5% oxygen.  
B - Not for use in atmospheres immediately dangerous to life or health.  
C - Do not exceed maximum use concentrations established by regulatory standards.  
J - Failure to proper use and maintain this product could result in injury or death.  
M - All approved respirators shall be selected, fitted, used, and maintained  
in accordance with MSHA, OSHA and other applicable regulations.  
N - Never substitute, modify, add, or omit parts. Use only exact replacement  
parts in the configuration as specified by the manufacturer.  
O - Refer to user's instructions, and/or maintenance manuals for  
information on use and maintenance of these respirators.  
P - NIOSH does not evaluate respirators for use as surgical masks.  
S - Users cannot wear the FFR without the silicone seal.

**SERVICE TIME LIMITATION**

1. Do not use the same respirator for more than 8 hours, continuous or intermittent,  
in dirty workplaces that could result in high filter loading.

2. Leave contaminated area and remove respirator if breathing becomes difficult or  
proper fit cannot be obtained or the respirator is damaged, soiled, or distorted.  
Do not use the same respirator again once it has been used in contaminated area.

3. The maximum service time of this respirator is five (5) consecutive calendar days  
(including days of non-use), beginning from the first day of use.

**LIMITED WARRANTY  
IMPORTANT NOTICE TO  
PURCHASER:**

This limited warranty is made in lieu of the warranties of  
merchantability, fitness for particular purposes and all  
other warranties, expressed or implied. There are no  
other warranties which extend beyond the description on  
the face hereof. The physical standards and specifica-  
tions of Makrite will be met by products sold.




Exclusive Remedies: damages for the breach of this  
limited warranty are limited to the replacement of such  
quantity of Makrite, products proved to be defectively  
manufactured. Except as provided above, Makrite shall  
not be liable or responsible for any loss, damage, or  
liability, direct, indirect, incidental, special, or consequen-  
tial, arising out of sale, use, or misuse, or the inability to  
use products by the user.

**WARNINGS:**

1. Wearer must read fitting instructions and use  
limitations before use.

2. This respirator is not guaranteed to protect against  
the risk of contracting a disease or infection.

3. Before occupational use of this respirator, a written  
respiratory protection program must be implement-  
ed meeting all the local government requirements.  
In the United States employers must comply with  
OSHA 29 CFR 1910.134 which includes medical  
evaluation, training, and fit testing.



Made in China

LOT NO.  
Expiration Date:



9800-N95  
TC-84A-XXXX  
20pcs/pack  
3 pcs Silicone Seal

Silicone Seal Respirator N95  
95% Efficient Protection Against  
Non-Oil Particulates



Respirator can not be used without the silicone seal



9800-N95  
TC-84A-XXXX  
20pcs/pack  
3 pcs Silicone Seal


Silicone Seal Respirator N95  
95% Efficient Protection Against  
Non-Oil Particulates




Respirator can not be used without the silicone seal

RESPIRATOR FITTING INSTRUCTION


- 1


- STEP 1**


After install the silicone seal (please refer to the Silicone Seal Installation Method below), hold the respirator in hand at your fingertips, allowing the headbands to hang freely below your hand.
- 2


- STEP 2**


Press the respirator firmly against your face on the bridge of your nose.
- 3


- STEP 3**

Stretch and position the top band high on the back of the head. Stretch the bottom band over the head and position below your ears.
- 4


- STEP 4**

Using both hands, to adjust the respirator to your nose.
- 5

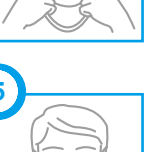
- STEP 5**

To test fit: a) Cup both hands over the respirator being careful not to disturb position, and b) inhale vigorously. If air leaks around the edges, tighten the upper headband by headband adjuster for better fit. Please carefully follow these fitting instructions during each use to achieve proper fit.

Silicone Seal Cleaning & Installation Method:

Cleaning: After each use, remove the silicone seal from the respirator and set respirator aside, and then to clean the silicone seal with water, wipe with 75% concentration alcohol or boil at high temperature .

These cleaning or disinfection methods will not affect the tightness and functionality of the product. It is recommended that the silicone seal can be reused 30 times.

High temperature boiling disinfection(only for the silicone seal): After boiling at 100 °C for 30 times, the silicone seal will not deform and affect the tightness. After cleaning or disinfection, place it in a ventilated and clean place, and install the silicone seal on the mask when it is completely dry.

- Installation: 1. Align the arrow symbol on the mask according to the direction indicated by the floating arrow on the silicone seal (as shown in Figure 1).
2. Put the silicone seal along the edge of the mask (as shown in Figure 2).
3. When the floating arrow on the silicone seal is aligned with the mask, it means that the mask is installed correctly.
4. When the silicone seal is combined, please ensure that the edge of the silicone rubber and the mask are smoothly and well attached.

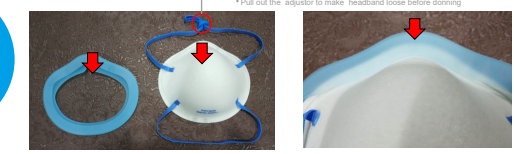


Figure 1

Figure 2

200 MM