

RESPIRATOR FITTING INSTRUCTION

1



STEP 1

After install the silicone seal (please refer to the Silicone Seal Installation Method below), hold the respirator in hand at your fingertips, allowing the headbands to hang freely below your hand.

2



STEP 2

Press the respirator firmly against your face on the bridge of your nose.

3



STEP 3

Stretch and position the top band high on the back of the head. Stretch the bottom band over the head and position below your ears.

4



STEP 4

Using both hands, to adjust the respirator to your nose.

5



STEP 5

To test fit: a) Cup both hands over the respirator being careful not to disturb position, and b) inhale vigorously. If air leaks around the edges, tighten the upper headband by headband adjustor for better fit.

Please carefully follow these fitting instructions during each use to achieve proper fit.

Silicone Seal Cleaning & Installation Method:

Cleaning: After each use, remove the silicone seal from the respirator and set respirator aside, and then to clean the silicone seal with water, wipe with 75% concentration alcohol or boil at high temperature.

These cleaning or disinfection methods will not affect the tightness and functionality of the product. It is recommended that the silicone seal can be reused 30 times.

High temperature boiling disinfection (only for the silicone seal): After boiling at 100 °C for 30 times, the silicone seal will not deform and affect the tightness. After cleaning or disinfection, place it in a ventilated and clean place, and install the silicone seal on the mask when it is completely dry.

Installation: 1. Align the arrow symbol on the mask according to the direction indicated by the floating arrow on the silicone seal (as shown in Figure 1)

2. Put the silicone seal along the edge of the mask (as shown in Figure 2)

3. When the floating arrow on the silicone seal is aligned with the mask, it means that the mask is installed correctly.

4. When the silicone seal is combined, please ensure that the edge of the silicone rubber and the mask are smoothly and well attached.

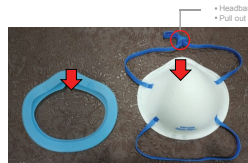


Figure 1



Figure 2

* Headband adjustor must be in-place to don the respirator and cannot be removed
* Pull out the adjustor to make headband loose before donning