


RESPIRATOR FITTING INSTRUCTION



STEP 1

Hold the respirator in hand with the nosepiece at your fingertips, allowing the headbands to hang freely below your hand.

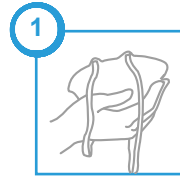


STEP 2

Press the respirator firmly against your face with the nosepiece on the bridge of your nose.



STEP 3

Stretch and position the top band high on the back of the head. Stretch the bottom band over the head and position below your ears.

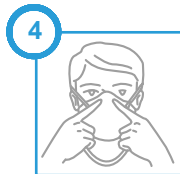

STEP 4

Using both hands, mold nosepiece to the shape of your nose.



STEP 5

To test fit: a) Cup both hands over the respirator being careful not to disturb position, and b) inhale vigorously. If air leaks around the edges, reposition the straps or adjust strap tension where applicable for better fit. Please carefully follow these fitting instructions during each use to achieve proper fit.

Caution:

The respirator should be fit checked prior to each wearing.



Headband adjuster must be in-place to don the respirator and cannot be removed
Pull out the adjuster to make headband loose before donning