

Instructions Manual

Filtering Facepiece Makrite

Please read this manual carefully before use. Keep it for reference.


Particulates Filtering Facepiece Makrite

is made of materials of polypropylene and meltbrown fabric.

As a private protection equipment man can use it against dust and particulates efficiently and comfortably.

Protection: N95-Particulate Filter (95% filter efficiency level). Effective against particulate aerosols free of oil (with diameter more than $\varnothing 0.3 \mu\text{m}$). Time use restrictions.

Application scope: extensive industry fields such as: building site, mining, spinning, grinding, pharmacy, cement, glass, hardware, casting...

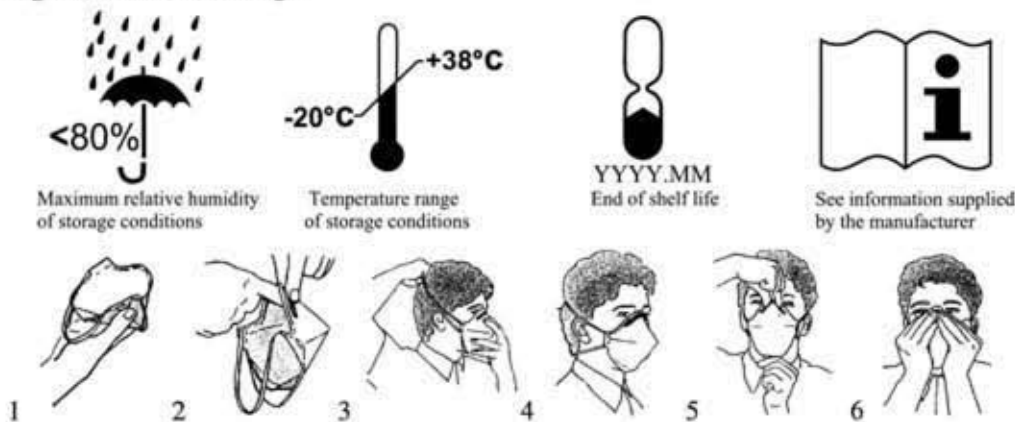

Direction for fitting:

Cup the mask in the hand, allowing the head-bands to hang freely.

Hold the mask under the chin with the nose-piece facing outwards.

Place the lower head-band around the neck below the ears. Holding the mask against the face with one hand, place the top head-band above the ears, around the crown of the head.

Mould the nose area to the shape of the face, running the fingertips of both hands from the top of the nose-piece down both sides while pressing inward. Check for leak-tightness (Facial fit) as follows:

- place both hands over the mask and inhale sharply
- negative pressure should be felt inside the mask
- if you feel air entering around the edges, re-adjust the mask by tightening the nose-piece and/or increasing tension of the head-bands
- repeat the procedure until satisfactory face-fit has been achieved.

Caution and limitation:

1. Not for use in atmospheres containing less than 19.5% oxygen
2. Not for use in atmospheres immediately dangerous to life or health.
3. Do not exceed maximum use concentrations established by regulatory standards.
4. Failure to properly use and maintain this product could result in injury or death.
5. All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA and other applicable regulations.
6. Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
7. Refer to users instructions, and/or maintenance manuals for information on use and maintenance of these respirators.
8. Always check the filtering facepiece before use. Be sure:
 - 1) No damage or deformation
 - 2) No holes in the breathing zone
 - 3) No mildewed points or unusual smellOtherwise the products cannot be used anymore.

Product Information:

Name: N95 filtering facepiece

Type: Makrite

Used material: Polypropylene and meltbrown fabric

Protection level: N95

NIOSH Approval Number:

Manufacturer: Makrite Industries Inc.