
MADE IN USA  
FABRIQUÉ AUX É.U.  
HECHO EN USA  
See Inside for CBRN Warning  
Single Use/Disposable  
8 Respirators  
Usage Unique/Jetable  
8 Respiradores  
Vea Adentro las Advertencias Sobre CBRN  
Desechable/Un Solo Uso  
8 Respiradores

4800  
PARTICULATE RESPIRATOR  
RESPIRADOR PARA PARTICULAS  
Mediano/Grande  
Moyen/Grand
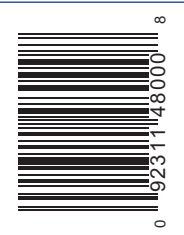
#### WARNING TO USER:

- Follow all instructions and warnings on the use of this respirator and wear during all times of exposure. Failure to do so will reduce respirator effectiveness, wear protection and may result in sickness or death.
- Before use, the user must first be trained by the employer in proper respirator use, in accordance with applicable safety and health standards, for the contaminant and exposure level in the assigned work area.
- The particulates and other contaminants which can be dangerous to your health include those which you cannot see, taste, or smell.
- This product has not been sold with warnings or use instructions for personnel involved in healthcare or related situations, where there may be the possibility of contact with disease or biological hazards. If you are considering such uses, first call the Moldex Technical Service Dept., +1 (310) 837-6500 or +1 (800) 421-0668, ext. 512/550.

#### RESTRICTIONS:

- This respirator does not supply oxygen and must not be used in atmospheres containing less than 19.5% oxygen.
- Do not remain in contaminated area if any physical distress occurs, for example breathing difficulty, dizziness or nausea.
- Leave contaminated area and replace respirator if it is damaged, distorted, a proper fit can not be obtained or breathing becomes difficult.
- Prior to each use, carefully inspect the entire respirator, including filter media and strap/weld attachment area for tears and damage. Staple or strap/weld attachment perforations do not affect respirator performance.
- Do not alter, modify, or abuse this respirator.
- Store unused respirators in box/bag in a clean, dry, non-contaminated area.
- Dispose of respirator according to your employer's policy and local regulations.
- Dispose no later than thirty days after first use.
- Use respirator before the "use by" expiration date printed on box/bag.
- If used for welding, wear appropriate eye and face protection.

For technical assistance call Moldex Technical Service Department, +1 (310) 837-6500 or +1 (800) 421-0668, ext. 512/550.

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#### LIMITED WARRANTY - IMPORTANT NOTICE TO PURCHASER


This limited warranty is made in lieu of the warranties of merchantability, fitness for particular purposes and all other warranties, express or implied. There are no other warranties which extend beyond the description on the face hereof. The physical standards and specifications of Moldex will be met by products sold. **Exclusive Remedies:** damages for the breach of this limited warranty are limited to the replacement of such quantity of Moldex products proved to be defectively manufactured. Except as provided above, Moldex shall not be liable or responsible for any loss, damage, or liability, direct, indirect, incidental, special, or consequential, arising out of sale, use, or misuse, or the inability to use products by the user.

#### USE AGAINST:

- Particulate aerosols free of oil.
- Nuisance levels (less than OSHA Permissible Exposure Limits) of Organic Vapor Odors and/or Ozone.

#### DO NOT USE AGAINST:

- Concentrations of contaminants which are unknown or are immediately dangerous to life or health.
- Concentrations of particulates which exceed the maximum use concentration or 10 times the OSHA Permissible Exposure Limit, whichever is lower.
- Gases, vapors, asbestos, paint spray, sandblasting or particulate materials which generate harmful vapors. This respirator is not approved for any vapors or gases, so when used against organic vapor odors and/or ozone do not exceed the OSHA Permissible Exposure Limit.
- Oil-based mists.
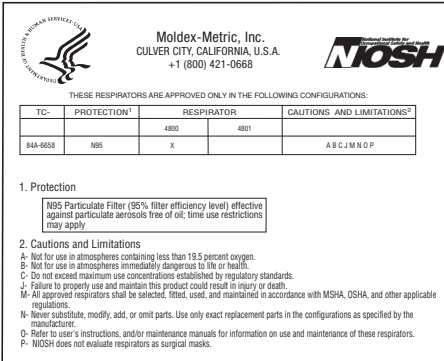
#### MOLDEX-METRIC, INC.

10111 W. Jefferson Blvd., Culver City, CA 90232  
Tel: +1 (800) 421-0668 or +1 (310) 837-6500  
Fax: +1 (310) 837-9563  
Email: sales@moldex.com www.moldex.com

#### AUSTRALIA - MOLDEX-METRIC, INC.

Unit 30/14-16 Stanton Road, Seven Hills, NSW 2147  
TEL +61 (02) 9638 7588 FAX +61 (02) 9638 7544  
Email: moldexaustralia@moldex.com www.moldex.com

#### NEW ZEALAND - MOLDEX-METRIC, INC.

3 Spence Road, Henderson, Auckland, NZ  
TEL +64 (0) 837 8247 FAX +64 (0) 837 8248  
Email: moldexnz@moldex.com www.moldex.com

Moldex Technical Service Department: +1 (800) 421-0668 or +1 (310) 837-6500 ext. 512/550

Moldex, Ideas that wear well, Dura-Mesh, Softspan and the PVC-Free logo are registered trademarks. The mesh design, waved design, AirWave, FlexFit are trademarks of Moldex-Metric, Inc. Patents pending.

Distributed by:

Printed on Recycled Paper

Moldex-Metric International LLC • AEN 64 964 843 487  
10111 W. Jefferson Blvd., Culver City, CA 90232, USA  
This warranty is given by  
Moldex will pay your reasonable direct expenses of making a claim under this limited warranty which is upheld. You may submit details  
take the product with proof of purchase to the store where you purchased the product to contact Moldex.  
Any claim under the warranty must be made within 1 year of the date of purchase of this product. To make a claim under the warranty,  
the ability to use products by the user, although you may have other rights and remedies under the Australian Consumer Law,  
not be liable or responsible for any loss, damage or liability, direct, indirect, incidental, special, or consequential, arising out of sale, use, or misuse, or  
physical standards and specifications of Moldex will be met by products sold. **Exclusive Remedies:** damages for the breach of this limited warranty are limited to  
replacement of the number of Moldex products proved to have been defectively manufactured. Under this limited warranty, Moldex shall  
to time, or other similar legislation of a state or territory of Australia (Australian Consumer Law). The product sold will meet Moldex's  
limited warranty or those in the Australian Consumer Law in Schedule 2 to the Competition and Consumer Act 2010 as amended from time  
accepted quality and the failure does not amount to a major failure. You are also entitled to have the goods repaired or replaced if the goods fail to be  
for other reasonably foreseeable loss or damage. You are also entitled to have the goods repaired or replaced if the goods fail to be  
This limited warranty is made in lieu of all other warranties, express or implied. There are no other warranties which extend beyond this  
cannot be excluded under the Australian Consumer Law. You are entitled to replacement or refund for a product failure and to compensation  
Our goods come with guarantees which cannot be excluded under the Australian Consumer Law. Our goods come with guarantees which cannot be excluded under the Australian Consumer Law.

#### LIMITED WARRANTY - IMPORTANT NOTICE TO AUSTRALIAN PURCHASER

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## 4800 Series

### N95 PARTICULATE RESPIRATOR

Plus Relief From Nuisance Levels (less than OSHA Permissible Exposure Limit) of Ozone and Organic Vapor Odors.

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## P2

Premium Welding Respirator  
Plus Relief from Nuisance Levels (Below Australian Exposure Standards) of Organic Vapours and Ozone.

- Exclusive Adjustable SmartStrap®
- Easy Breathing Waves
- Ultimate FlexFit® Comfort

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
8 Respirators  
8 Respirateurs  
8 Respiradores

#### FITTING

- Users must follow instructions each time respirator is worn.
- It's the wearer's responsibility to obtain a proper fit.
- To choose the correct model/size Moldex Respirator you must fit test and user seal check.
- If you cannot obtain a proper fit do not enter the contaminated area.
- OSHA regulation 1910.134 (f) requires that the user must be fit tested before being assigned a respirator and annually thereafter.
- Use Moldex Qualitative Fit Test Kit. User must be clean-shaven. Any facial hair, such as beards or long sideburns, may prevent the respirator from fitting properly.
- Eye wear must not interfere with face to facepiece seal.
- Re-fit test when there are any changes to the wearers physical condition that may affect fit.

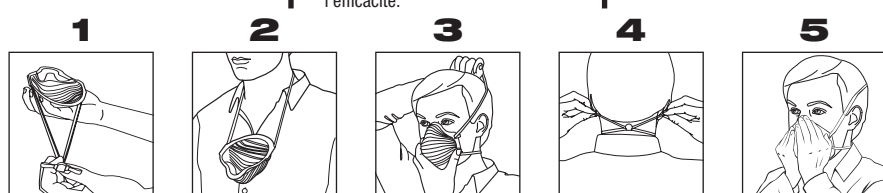- Pull adjustment clip and strap fully below the bottom of the mask. (Fig 1)
- Place strap around neck so outside of shell is against chest. (Fig 2)
- a) Fit mask to face and pull top of strap to crown of head. Nose cushion must not be creased inside of respirator.  
b) Always use both top and bottom straps. (Fig 3)
- To tighten pull both strap tabs. (Fig 4)
- To loosen, grasp both sides of neck strap and pull.
- Each time user enters contaminated area respirator must be seal checked. Cover front of respirator by cupping both hands. (Fig 5) INHALE SHARPLY. A negative pressure should be felt inside respirator. If any leakage is detected at respirator edges, tighten or adjust strap by pulling back along the sides and/or reposition respirator. Repeat until sealed properly, otherwise see your supervisor. Entry into a contaminated area with an improper fit may result in sickness or death.
- To hang around neck, pull top strap away from head, then pull mask down. (Fig 2)

**CAUTION:** When changing from any model/size to another respirator, you must fit test. If used against biohazards discard immediately after removing from face and do not hang around neck.

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Failure to follow all instructions and warnings on the proper use of this respirator may reduce respirator effectiveness, wear protection and may result in sickness or death. For additional information consult your supervisor, read all information on box, and phone +1 (800) 421-0668 or +1 (310) 837-6500, ext. 512/550.

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