


# Instruction for Use:

## To Put On This Respirator:

1. Hold the respirator in your hand with the metal adjusting strip away from you, allowing the headbands to hang below your hands.
2. Press the respirator firmly against your face with the narrow end on the bridge of your nose.
3. Pull the top strap over your head and position it high on the back of your head. Pull the bottom strap over your head and position it below your ears.
4. Using both hands, mold the metal strip to the shape of your nose.
5. Check the facepiece-to-face seal by cupping both hands over the respirator and exhale vigorously. If air escapes around your nose, reform the metal strip. If air escapes around the edges, reposition the straps. Repeat as necessary until no air escapes, and you have obtained a satisfactory seal.

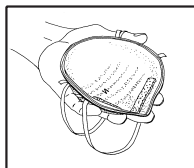

## Time Use Limitations:

The use and reuse of this respirator is subject to considerations of hygiene, damage and increased breathing resistance. For dirty workplaces, service time should not be extended beyond 8 total hours of use.