


Instructions for fitting – Model F550G (Duckbill shape)

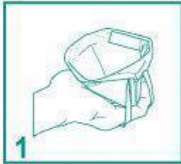

1. Prepare to use respirator

- Unfold the respirator by pulling top part (with foam) up and bottom part down so the respirator is opened all the way.
- Check the respirator for holes, tears or other damage to surfaces, straps, and nose foam to ensure that your respirator works well. Replace with a new respirator if damaged.
- Place the mask within in hand, the nosepiece pointing at your fingertip and the headbands hanging downwards.

If you wear glasses, read steps 2 through 4 then remove your glasses before putting on the respirator.



2. Place respirator on face

- Place the respirator against your face, with the bottom under your chin and nose foam across bridge of your nose.
- Pull the top elastic band over your head and place high at the top back of your head.
- Next, pull the bottom elastic band over your head and place around the neck below your ears.



3. Adjust the nosepiece

- Use both hands to bend the nosepiece to fit snugly against your nose and face. The respirator may not fit as well if you pinch the nosepiece using one hand. Use two hands.
- Slide fingers down both sides of nosepiece to seal it against your nose and face.



4. Check the face seal

- Completely cover the middle panel of the respirator with one or both hands.
- Do not push the respirator against your face.
- With your hands in place on the surface of the respirator, exhale, or breathe out sharply.
- If you feel air blowing on your face or eyes, the respirator needs to be adjusted. To adjust, stretch the headbands again by pulling them backwards. Then correct the position of the respirator and check the face seal again. When the respirator is a good fit, you will not feel any air blowing on your face or eyes.

Return eyewear to face, if applicable.

Proper fit of the respirator is critical to reducing wearer exposure to airborne germs. It is important to read User Instructions carefully.

S-Special or Critical User's Instruction

The model PASTURE F550G filtering facepiece respirator has been manufactured by Approval Holder Champak Enterprise Co., Ltd. for Pasture Pharma PL under approval TC-84A- XXXX

Qty: 20 pcs