

**WARNING**


This respirator helps protect against certain particles. Misuse may result in sickness or death. For more information please call priMED at 1-877-877-4633.

- IMPORTANT:
- Before use, the wearer must read and understand these user instructions
  - Contains no natural rubber latex
  - NIOSH does not evaluate respirators as surgical masks
  - Before use, the wearer must be fit tested according to 29CFR 1910.134

NIOSH APPROVED N95 RESPIRATOR	USE FOR:	DO NOT USE FOR:
	<ul style="list-style-type: none"><li>• Solids such as those from processing minerals, and other substances.</li><li>• Liquids or non oil-based particles from sprays that do not also emit oil aerosols or vapors.</li><li>• Non-visible particles that may pose a health hazard</li></ul>	<ul style="list-style-type: none"><li>• Paint spray, oil aerosols, gases, vapors, asbestos, or sandblasting.</li></ul>
95% efficient against particulate aerosols free of oil		

CAUTIONS	LIMITATIONS
<ul style="list-style-type: none"><li>• Failure to follow all instructions may reduce respirator effectiveness.</li><li>• Leave the contaminated area immediately if dizziness, irritation or other distress occurs.</li><li>• Store the respirator in a clean, dry environment, away from contaminated areas.</li><li>• Discard used respirator in accordance with applicable regulations.</li></ul>	<ul style="list-style-type: none"><li>• This respirator does not supply oxygen. Do not use in atmospheres containing less than 19.5% oxygen.</li><li>• Do not use where concentrations of contaminants are immediately dangerous to life and health.</li><li>• Do not alter or misuse the respirator.</li><li>• Do not use with beards or other facial hair or conditions that prevent a good seal between the face and the edge of the respirator.</li></ul>

**TIME USE LIMITATIONS:**  
If respirator becomes damaged, soiled, or breathing becomes difficult, leave the contaminated area and replace the respirator.

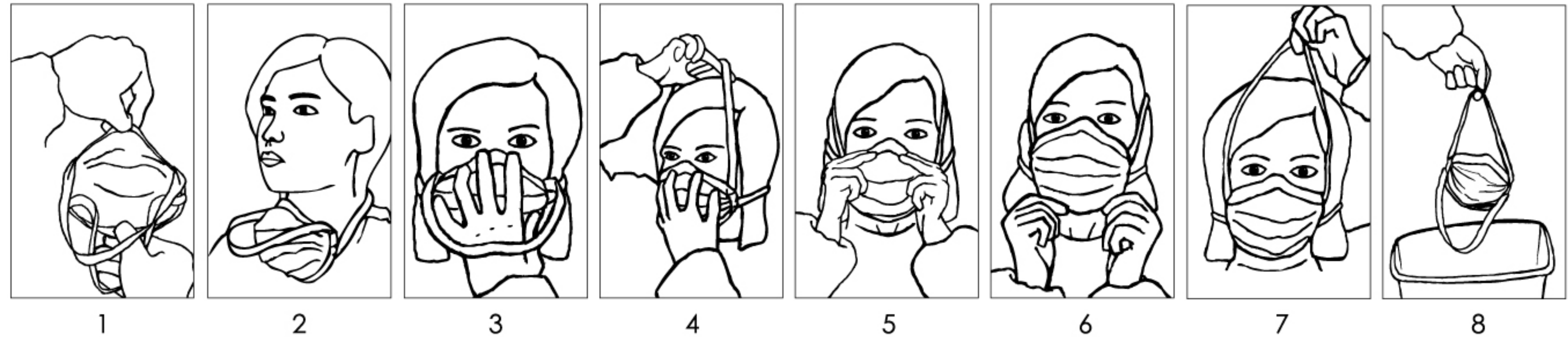

- INSTRUCTIONS FOR USE:**
- 1) Open the pleats. Bend the nose piece into a curved shape, centered in the middle.

2) Pull bottom strap over head to below ears.

3) Hold the mask in position over the nose and mouth with one hand.

4) With the other hand pull the top strap up to rest on the crown of your head. Ensure the straps are not twisted, but lie flat.

5) Using both hands, form the nose piece around your nose, and down over your cheek bones.

6) Perform a seal check: cup both hands around the perimeter of the mask then inhale and exhale sharply. If any air flows around the nose, tighten the nose piece. If air flows around the edges, reposition and adjust the piping.

7) Remove by the strap or edges.

8) Discard into trash bin.



priMED Medical Products, Inc.  
Edmonton, Alberta Canada  
1-877-877-4633



**THIS RESPIRATOR IS APPROVED ONLY IN THE FOLLOWING CONFIGURATION:**

TC-	Protection <sup>1</sup>	Respirator			Cautions and Limitations <sup>2</sup>
		PN4-7064	PN4-7074	PN4-7084	
84A-PMM4a	N95	X	X	X	ABCJMNOP

- 1. Protection**
- N95 - Particulate Filter (95%) filter efficiency level)  
Effective against particulate aerosols free of oil; time use restrictions may apply

- 2. Cautions and Limitations**
- A - Not for use in atmospheres containing less than 19.5% oxygen.
  - B - Not for use in atmospheres immediately dangerous to life or health.
  - C - Do not exceed maximum use concentrations established by regulatory standards.
  - J - Failure to properly use and maintain this product could result in injury or death.
  - M - All approved respirators shall be selected, fitted, used and maintained in accordance with MSHA, OSHA and other applicable regulations.
  - N - Never substitute, modify, add or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
  - O - Refer to users instructions, and/or maintenance manuals for information on use and maintenance of these respirators.
  - P - NIOSH does not evaluate respirators for use as surgical masks.

PART# PN4-7064, PN4-7074, PN4-7084



priMED

MANUFACTURE • DELIVER • PROTECT

priMED Medical Products Inc.

3<sup>rd</sup> Floor, 1259-91 Street SW

Edmonton, AB, Canada T6X 1E9

www.primed.ca 1-877-877-4633

September 2011