



DS-SM-P2W



Shanghai Dasheng Health Products Manufacture Co., Ltd

Phone : +86-21-57783126 Website: www.dashengmask.com

FITTING INSTRUCTIONS:

To be followed each time respirator is worn.



1. Hold the preformed mask in the palm of your hand, with the nose bar towards the top and the straps hanging down.



2. Place the mask under your chin, covering your mouth and nose.



3. Take the fixing device on which the elastic straps are linked, pass it around the back of your head. Fasten the fixing device through snap.



4. Adjust the mask onto your face




5. Perform a User Seal Check prior to each wearing. To check fit, close off the inlet opening of the exhalation valve by covering with the palm of the hand, inhale gently so that the face piece collapses slightly, and hold the breath for ten seconds. If the face piece remains in its slightly collapsed condition and no inward leakage of air is detected, the tightness of the respirator is considered satisfactory. If you CANNOT achieve a proper fit, DO NOT enter the contaminated area. See your supervisor.

REMOVAL INSTRUCTIONS:

See step 3 of Fitting Instructions and cup respirator in hand to maintain position on face. Pull bottom strap overhead. Still holding respirator in position, pull top strap over head and remove respirator.

This respirator contains no components made from natural rubber latex.

REPLACE FILTER INSTRUCTIONS:

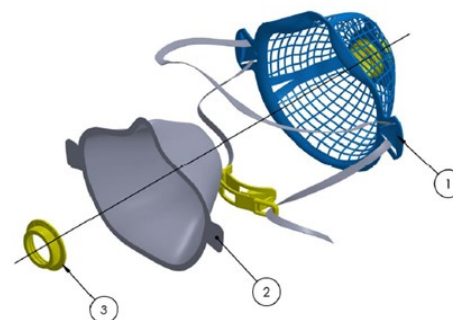

Steps to replace the replaceable Filter ②:

- 1 – Unscrew the Reusable Filter Lock ③ from the Reusable Filter Frame ①
- 2 – Take off the replaceable Filter ②
- 3 – Put a new replaceable Filter ② into position
- 4 – Screw back the Reusable Filter Lock ③ onto the Reusable Filter Frame ①



Shanghai Dasheng Health Products Manufacture Co., Ltd
No.228 Shihui Road, Songjiang, Shanghai 201613, China
Tel: +86-21-57783126 Fax: +86-21-57784148



These respirators are approved only in the following configurations:

Respirator Components					
TC-	Protection ¹	Filter Frame	Filter Lock	Filter	Cautions and Limitations ²
		RFF	RFL	DS-SM-FFP2W	
84A-8635	N95	X	X	X	ABCJMNOPS

1. PROTECTION:

N95 - Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil; Time use restrictions may apply.

2. CAUTIONS AND LIMITATIONS:

- A. Not for use in atmospheres containing less than 19.5 percent oxygen.
- B. Not for use in atmospheres immediately dangerous to life or death.
- C. Do not exceed maximum use concentrations established by regulatory standards.
- J. Failures to properly uses and maintain this product could result in injury or death.
- M. All approved respirators shall be selected, fitted, used and maintained in accordance with MSHA, OSHA and other applicable regulations.
- N. Never substitute, modify, add or omit parts. Use only exact replacement parts in the configuration as specified by the configuration as specified by the manufacturer.
- O. Refer to user's instructions and or maintenance manuals for information on use and maintenance of these respirators.
- P. NIOSH does not evaluate respirators as surgical masks.
- S. Special or critical User Instructions and/or specific use limitations apply. Refer to User Instructions before donning.

NIOSH APPROVED: N95

95% efficient against particulate aerosols free of oil. This respirator helps protect against certain particulate contaminants but does not eliminate exposure to or the risk of contacting disease or infection. Misuse may result in sickness or death.

S - Special or Critical User's Instructions: This respirator offers nuisance level relief from organic vapors not exceeding the US OSHA PEL or applicable government occupational exposure limits, whichever is lower.

FOR MORE INFORMATION and assistance on

Shanghai Dasheng Health Products Manufacture Co., Ltd

Phone: +86-21-57783126 or Website: www.dashengmask.com