




# SPECTRA SHIELD™ PROTECTOR N95

## *Effective Respiratory Protection*


Users shall follow the below fitting instructions each time the respirator is worn. This respirator is designed to be used as part of a respiratory protection program. This respirator must be used in accordance with OSHA Regulations 29 CFR 1910.134 and specifically to the OSHA requirement for fit testing before use. Users of this respirator must be clean-shaven and free of any facial hair, which may interfere with seal area and fit of this respirator. The users shall be fit tested with the same make, model, style, and size respirator as to be used.

### Step 1

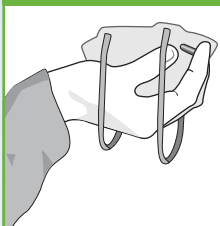


Hold the respirator in hand with nosepiece at your fingertips, allowing the head bands to hang freely below your hand.

### Step 2

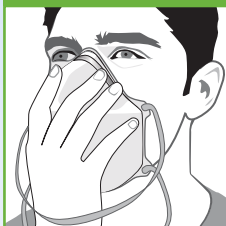

Place under chin and press the respirator firmly against your face with the nosepiece on the bridge of your nose.

### Step 3

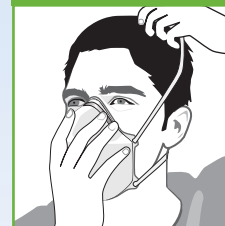


Stretch and position the top band high on the back of your head. Stretch the bottom band over the head and position it below your ears.

### Step 4



Using both hands, mold the metal nosepiece to the shape of your nose.

### Step 5

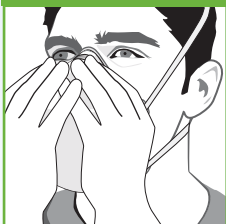

To test fit, cup both hands over the respirator and exhale vigorously. If air flows around your nose, tighten the nosepiece. If air leaks around the edges of the respirator, reposition the straps.

### Step 6



Repeat adjustments until the respirator is sealed properly. Entry into a contaminated area with an improper fitting respirator may result in illness or death.

Made in the U.S.A.  
Assembled in China.

**Special or Critical User's Instructions:** The model SpectraShield™ Protector N95 filtering facepiece respirator has been manufactured by approval holder Foss Manufacturing Company, LLC for Nexera Medical, Inc. under TC-84A-5670.

## REMEMBER:

Careful observance of these fitting instructions is an important step in safe respirator use.