## INSTRUCTIONS / INSTRUCCIONES

**To Put On This Respirator: Mise en place du respirateur: Para ponerse este respirador:****"S"- Special Instructions Section**

The RPD725N95 Air Purifying Filtering Facepiece N95 respirator has been manufactured by Venus Safety & Health Pvt. Ltd. For DSI Safety Inc. under TC-84A-8126. NOTE: Before issuing any respirator to be worn, a face fit test must be performed either qualitatively or quantitatively. See OSHA regulation 29 CFR 1910.134 or CSA Z94.4. Mask should not be worn for more than 8 hours.

**Pre-use instructions:**

- Before occupational use of this respirator, a written respiratory protection program must be implemented meeting all the local government requirements. In the United States, employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training and fit testing.
- Each time before putting on this respirator, inspect the facepiece to assure that there is no holes in the breathing zone.

**To Put On This Respirator:**

- Open the respirator by molding the top nose strip to the shape of your nose.
- With your chin in the bottom of the respirator, pull the bottom strap over your head and below your ears.
- Place the top of the respirator over your nose and pull the top strap over your head. Position the strap high on the back of your head. If straps are not tight enough to t well and offer a proper seal, slide the knob on the strap to shorten the straps, this will allow to tight the respirator on your face.
- Using both hands, mold the metal strip to the shape of your nose and position the respirator rmly against your face.
- Check the facepiece-to-face seal by cupping both hands over the respirator and inhale vigorously. A negative pressure should be detected inside the facepiece. If air leaks in from around the edges, reposition the straps. Repeat as necessary until no air leak is detected, and you have a satisfactory seal. If you cannot achieve a good face-to-facepiece seal do not enter the contaminated area and see your supervisor.

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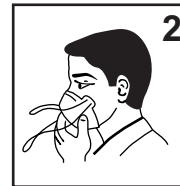
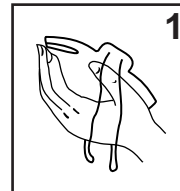
**Time Use Limitations:**

The use and reuse of this respirator is subject to considerations of hygiene, damage and increased breathing resistance. Replace if the respirator becomes damaged or if breathing becomes difficult.

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