
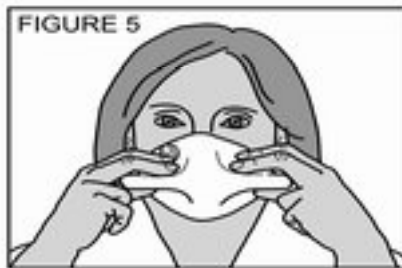
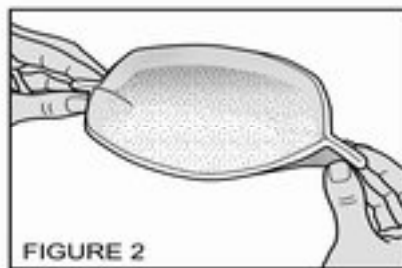


## Viramask™ N99ESC Model N99VSO7

**USER INSTRUCTIONS:** (Must be followed each time the VIRAMASKĬ is worn. **Do not** apply over full facial beards, cuts, abrasions or rashes.)

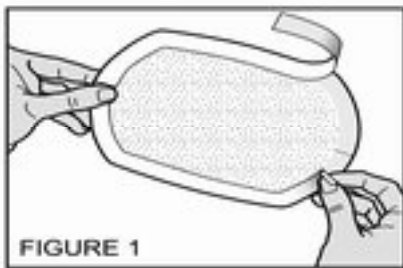

**Helpful Hints:** Applying a No Sting Barrier Film such as 3M Cavilon #3344, or any other barrier film to the face prior to applying the VIRAMASKĬ can help protect against damage or irritation to the skin during its removal. On removal, wrinkle nose and facial muscles while slowly removing the mask to break the adhesion bond quickly and easily.

**ATTENTION VIRAMASK™ USERS:** In accordance with OSHA requirements, the user should have the VIRAMASKĬ fit tested either quantitatively with a Portacount Plus Fit Tester or qualitatively.

### **APPLYING THE VIRAMASK™**

**Fig 1:** Peel off **ALL** of the protective adhesive release liner from the back of the mask exposing all of the adhesive.

**Fig 2:** Pre-pinch both sides of the mask, starting at the release liner slits, to align the adhesive.

**Fig 3:** Position the top of the mask across the bridge of the nose, just under the eyes, as shown and **press firmly** in place (for about 30 seconds).

**Fig 4:** Open mouth slightly to allow for talking, and position the bottom of the mask to the bottom of the chin, as shown and **press firmly** (for about 30 seconds).

**Fig 5:** With top and bottom of the mask in place, press the pinched portions firmly against the facial skin to provide a leak free seal in this area.

Continue pressing firmly ensuring there are no gaps or wrinkles so that the mask forms a smooth continuous surface to the skin around the perimeter of the mask. Following these instructions carefully creates a personalized custom face-seal fit to your particular facial structure while still allowing for ease of breathing and talking.

**User Seal Check:** According to OSHA manufacturer's allowable preference, with the hands covering most of the mask, positive pressure will displace the mask slightly forward away from the face when the user breathes out forcefully. If leaks are detected, apply more pressure to the area of the leak to more firmly set the adhesive.

### **REMOVING THE VIRAMASK™**

**Step 1:** With two fingers, press the skin next to the pinched section of the mask to prepare for its removal.

**Step 2:** With the other hand, **STARTING AT THE CHIN, PEEL UPWARD**, ōlow and slowō keeping it close to the skin surface to avoid skin damage. As the mask is removed; continue moving fingers as necessary, supporting the newly exposed skin by pressing on it while removing the mask.