




JACKSON^{*}
SAFETY Brand

64240

R10 Particulate Respirator with Valve N95

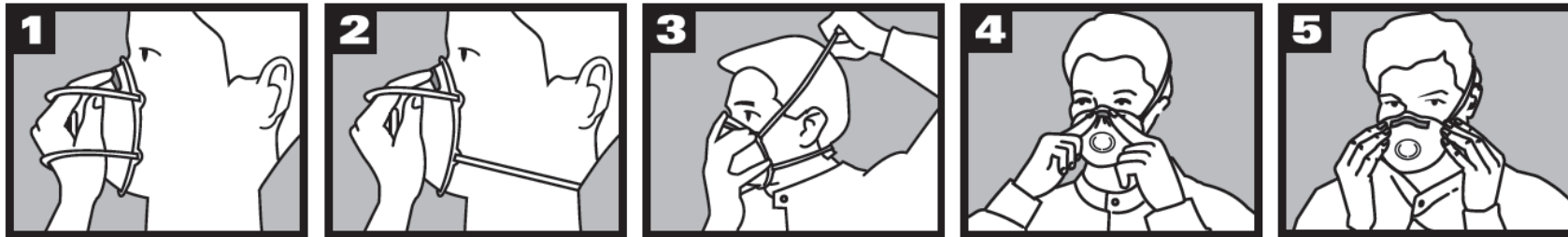

EN INSTRUCTIONS FOR USE:

Follow these instructions each time the respirator is worn. (See pictures 1-5)

1. Hold the respirator in your hand with the head straps hanging forwards over the front of the nose clip. Place the respirator over your nose and mouth. 2. Place the bottom elastic strap over your head just below the ears. 3. Place the top head strap over your head above your ears. 4. Using both hands, carefully shape the nose clip to fit your nose and adjust the respirator to achieve a completely tight fit to your face. 5. Test the tightness of the facial fit by covering the front of the respirator with both hands and inhale sharply. If air leaks out around the nasal area (nose clip) or elsewhere, these leaks must be eliminated by readjusting the nose clip and respirator, before entering the contaminated area.