


## User's Instructions for San Hwei Model No. SH1850A

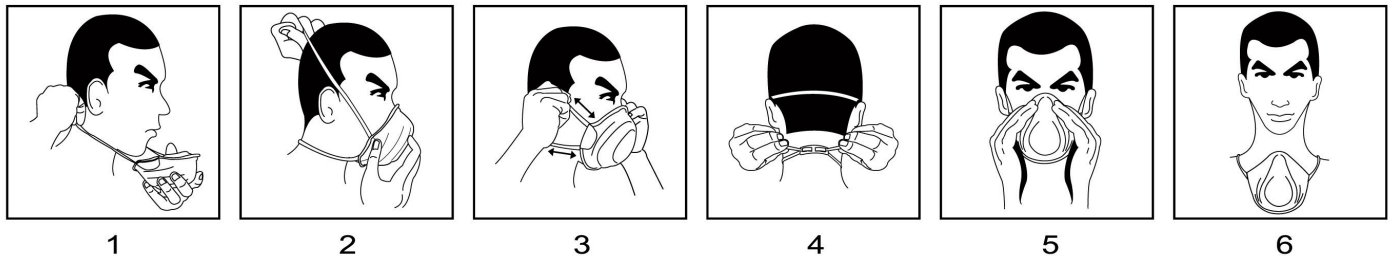

### **FITTING INSTRUCTIONS :**

1. Hold mask in position over nose and mouth, with retaining straps around back of head. Pull lower portion of strap over head and fit around neck.
2. Stretch and pull the upper strap over your head on the back of your head above your ears.
3. Adjust the tension in the upper and lower portions of the continuous loop strap, to ensure a positive face seal.
4. Tension may be increased by pulling strap at either side of buckle to obtain a tighter facial fit. Tension may be decreased by pushing out the back of the buckle to obtain a looser facial fit.
5. Cup both hands over the respirator and exhale vigorously. If air flows around your nose, reposition the respirator. If air leaks around the edges, repeat point 3 and 4.
6. During work break, take off the lower portion of strap. Let the respirator hang around the neck.
7. Change respirator immediately if breathing becomes difficult or respirator becomes damaged or distorted, or a proper face fit can not be maintained. Careful observance of these instructions is an important step in safe respirator use.

### **WARNING :**

1. This product does not eliminate the risk of contracting any disease or infection.
2. Failure to follow all instructions and limitations on the use of this product could reduce the effectiveness of respirator and result in sickness or death.
3. Before occupational use of this respirator, a written respiratory protection program must be implemented meeting all the local government requirements. In the United States, employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing.
4. This product does not supply oxygen. Use only in adequately ventilated areas containing sufficient oxygen to support life. Do not use this respirator when oxygen concentration is less than 19.5%.
5. Do not use when concentrations of contaminants are immediately dangerous to health or life.
6. Leave work area immediately and return to fresh air if (a) breathing becomes difficult, or (b) dizziness or other distress occurs.
7. Facial hairs or beards and certain facial characteristics may reduce the effectiveness of this respirator.
8. Never alter or modify this respirator in any way.
9. Only for single use. No maintenance necessary. Discard the used respirator after single use.
10. Keep respirators in the display box away from direct sunlight until use.
11. This respirator offers nuisance level relief from organic vapors (OV)/acid gas (AG) that are below the Permissible Exposure Limit (PEL). Nuisance level refers to concentrations not exceeding the OSHA PEL or other government occupational exposure limit, whichever is lower.

Rev. : A      Date : Oct. 2<sup>nd</sup>, 2017

*Note : The NIOSH approval label and user's instruction shall be supplied with the smallest commercial packaging. If the label or user's instruction is not printed on the outside of the container, it will be included as an insert inside the container.*