



## DONNING INSTRUCTIONS

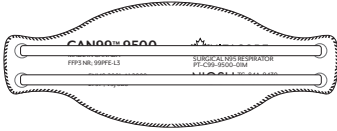


1

Remove respirator from its packaging positioning both straps in the front of the mask.

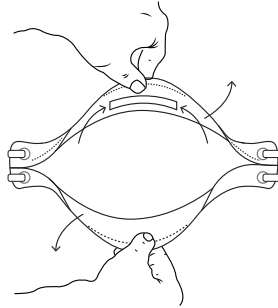

2

Fully open the top and bottom panels, gently bending the nose piece to form a round arch.

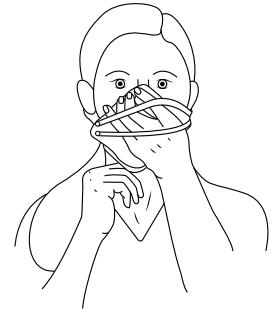


3

With one hand under the straps in front, place the respirator on your nose bridge with the bottom panel securely under your chin.



4

Pull the bottom strap over your head and place it under your neck below your ears. Pull the top strap over your head positioned high on the back of your head.

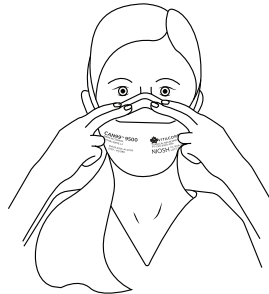

5

Before sealing nose bridge, ensure top panel is resting high on your nose bridge and bottom panel under your chin. Place fingertips from both hands on the top of the nose bridge and press down along the nosepiece molding it to the shape of your nose.



6


**PERFORM USER SEAL CHECK**  
Check the seal by placing both hands over the respirator and inhale and exhale sharply. If leakage is detected, then lift the nose bridge higher and repeat step 6. If you cannot achieve a proper seal do not use respirator in a contaminated area.

## DOFFING INSTRUCTIONS




1

When removing the respirator DO NOT TOUCH the facepiece. Gently lift the bottom strap around your neck over your head.

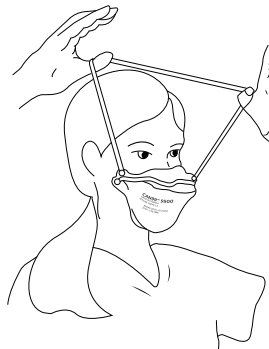

2

Lift the top strap without touching the facepiece.

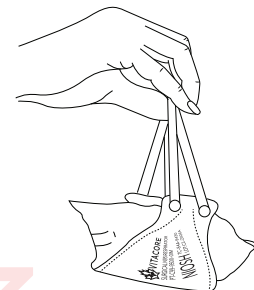

3

Dispose or recycle your respirator in accordance with your local infection control policies.

**Warning:** Respirators help protect against certain airborne contaminants. Before use, the wearer must read and understand the User Instructions provided as a part of the product packaging. Follow all local regulations. In the U.S., a written respiratory protection program must be implemented meeting all the requirements of OSHA 29 CFR 1910.134 including training, fit testing and jurisdiction, as appropriate. Misuse may result in sickness or death. For correct use, consult supervisor or User Instructions.